NGene**AnalySys™**

Instructions for Use

NGeneAnalySys™

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The screenshot shows the 'Pipeline Information' window for version v1.2, released on 2018-08-28. The 'Annotation Database' tab is selected. It displays a table of databases used in the pipeline:

Category	Database	Version	Release	Description	Source
Disease Database	ClinVar	201808	Jul.2018	Database of assertions about the clinical significance and phenotype relationship of human variations.	https://www.ncbi.nlm.nih.gov/clinvar
Disease Database (Somatic)	Catalogue of Somatic Mutations In Cancer (COSMIC)	v83	Dec.2017	Database for exploring the impact of somatic mutations in human cancer.	https://cancer.sanger.ac.uk/cosmic
Disease Database (Somatic)	Precision Oncology Knowledge Base (OncoKB)	v1.12	Oct.2017	Precision oncology knowledge base and contains information about the effects and treatment implications of specific cancer gene alterations.	http://oncokb.org
Human Reference	Mills	hg19	Dec.2013	An initial map of insertion and deletion (INDEL) variation in the human genome.	https://software.broadinstitute.org/gatk/down
Human Reference	Reference Genome Sequence	hg19	Feb.2009	Source of full human genome reference sequences.	https://www.ncbi.nlm.nih.gov/grc/human
Human Reference	RefSeqGene	.	Oct.2018	RefSeqGene defines genomic sequences to be used as reference standards for well-characterized genes and is part of the Locus	https://www.ncbi.nlm.nih.gov/refseq/rsg

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The screenshot shows the 'Pipeline Information' window for version v1.2, released on 2018-08-28. The 'Tools' tab is selected. It displays a table of analysis tools used in the pipeline:

Software	License	Version	Release	Description	Source
BBDuk	the University of California	v37.23	May.2017	Filters or trims reads for adaptors and contaminants using k-mers	https://github.com/lh3/bbtools/
BCFtools	MIT	v1.3.1	Apr.2016	Discovers short variants from Next Generation Sequencing data	http://www.htslib.org
BWA	GPLv3, MIT	v0.7.10	Jul.2014	Mapping low-divergent sequences against a large reference genome	http://bio-bwa.sourceforge.net
FastQC	GPLv3	v0.11.3	May.2015	Quality control tool for high throughput sequence data	http://www.bioinformatics.bbsrc.ac.uk/projects/fastqc
FastQValidator	the University of Michigan	v0.1.1	Oct.2012	Validation the format of FASTQ file	http://gerome.sph.umich.edu/wiki/FastQValidator
GATK	BSD 3-clause	v4.0.6.0	Jul.2018	Wide variety of tools with a primary focus on variant discovery and genotyping	https://software.broadinstitute.org/gatk
Pindel	GPL	v0.2.5b9	Jul.2016	Detect breakpoints of large deletions, medium sized insertions, inversions, tandem duplications and other structural variants at single-based resolution from next-gen sequence data	http://gmt.genome.wustl.edu/packages/pindel
SAMtools	BSD, MIT	v1.3.1	Aug.2016	Various utilities for manipulating alignments in the SAM format	http://www.htslib.org

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The screenshot shows the 'New Analysis Request' page. It includes fields for 'Run' name, 'Instrument' (Illumina MiSeqDX selected), 'FASTQ File Source' (Local Fastq Files selected), and detailed sample information for two entries (NA12878-2 and NA12878-1) across four categories: Sample, Panel, Source, and Disease. The 'Sample' column lists sample IDs, the 'Panel' column lists BRCAaccuTest PLUS, the 'Source' column lists Peripheral Blood, and the 'Disease' column lists HBOC. At the bottom are 'CANCEL' and 'SUBMIT' buttons.

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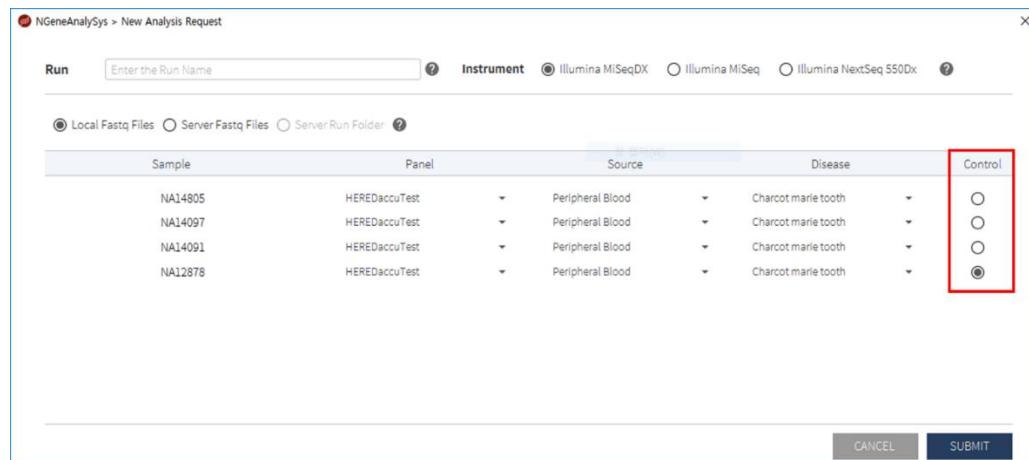
The screenshot shows the 'New Analysis Request' dialog box. At the top, there are tabs for 'Run' (selected), 'Instrument' (Illumina MiSeqDX), and other sequencing platforms. Below these are buttons for 'Local Fastq Files', 'Server Fastq Files', and 'Server Run Folder'. The main area is a table with columns: Sample, Panel, Source, Disease, and Control. The 'Control' column for the last row (NA12878) contains a radio button that is checked and highlighted with a red box. The table data is as follows:


Sample	Panel	Source	Disease	Control
NA14805	HEREDaccuTest	Peripheral Blood	Charcot marie tooth	<input type="radio"/>
NA14097	HEREDaccuTest	Peripheral Blood	Charcot marie tooth	<input type="radio"/>
NA14091	HEREDaccuTest	Peripheral Blood	Charcot marie tooth	<input type="radio"/>
NA12878	HEREDaccuTest	Peripheral Blood	Charcot marie tooth	<input checked="" type="radio"/>

At the bottom right of the dialog are 'CANCEL' and 'SUBMIT' buttons.

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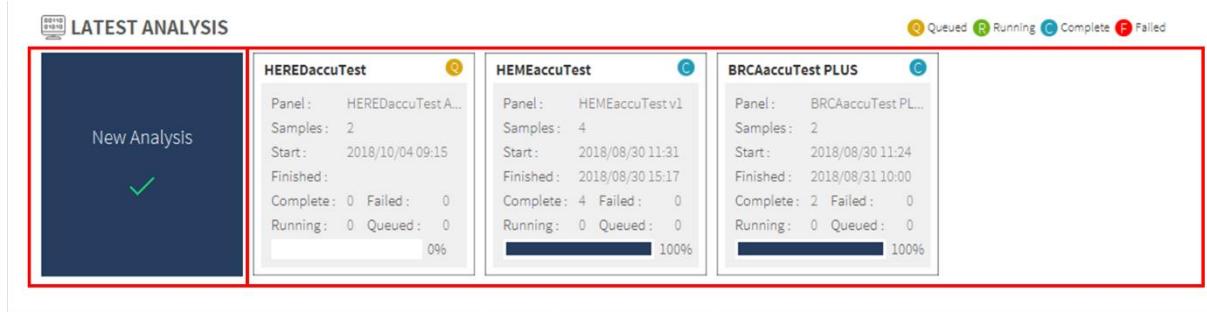Q Queued
 R Running
 C Complete
 F Failed

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The screenshot shows the NGeneAnalySys software interface. The top navigation bar has three tabs: DASHBOARD, RESULTS (which is highlighted with a red box and a yellow circle labeled '1'), and SAMPLES. Below the navigation bar are two tables of analysis results. The first table is for the HEMEaccuTest run, which used an Illumina MiSeqDx sequencer. It lists four samples: 15813AL-core-panel, Blood-22, Blood-38, and Blood-65. The second table is for the BRCAaccuTest PLUS run, also using an Illumina MiSeqDx sequencer. It lists two samples: NA13715-1 and NA14623-1. Both tables show columns for Run, Sequencer, Start Date, Finished Date, Sample, Panel, Status, Variant, and QC. A legend at the bottom indicates color coding for variant types: C (green), P (red), LP (orange), US (blue), LB (purple), and B (dark blue). On the right side, there is a search panel with fields for RUN, a search button, and a RESET button. A yellow circle labeled '2' points to the 'SEARCH' button. A yellow circle labeled '3' points to the 'RESULTS' tab. A yellow circle labeled '4' points to the 'EXCEL' icon in the search panel.

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The screenshot shows the NGeneAnalySys software interface. On the left, there's a table titled 'RESULTS' displaying analysis data for four samples: 15813AL-core-panel, Blood-22, Blood-38, and Blood-65. The columns include Run, Sequencer (Illumina MiSeqDx), Start Date (2018/08/30 11:31:23), Finished Date (2018/08/30 15:17:31), Sample, Panel (HEMEaccuTest v1), Status (C), Variant counts (T1, T2, T3, T4), and QC (WARNING, PASS). On the right, a search panel is open with the following fields:

- DATE**: Two date input fields labeled '2' and '3'.
- PANEL**: A dropdown menu set to 'HEMEaccuTest v1'.
- DATE**: A date range selector from '2018-01-01 ~ 2018-12-31'.

Buttons for **RESET** and **SEARCH** are at the bottom of the search panel. The status bar at the bottom right indicates 'System version : 1.4.6.0'.

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
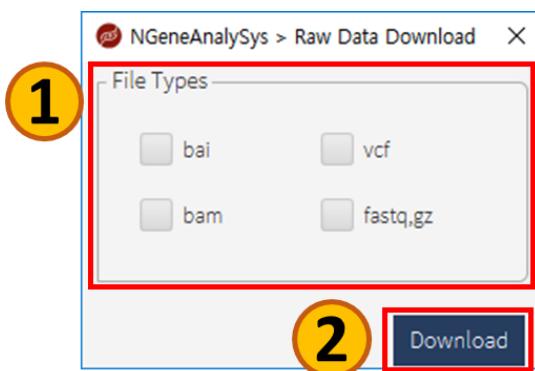
The screenshot shows the NGeneAnalySys software interface. The top navigation bar includes the logo, “NGeneAnalySys™”, “DASHBOARD”, “RESULTS” (which has a red box around it and a yellow circle with the number 1 above it), “SAMPLES” (which has a red box around it and a yellow circle with the number 2 above it), and “Test” dropdown. Below the navigation bar, the sample information is displayed: NA12878, RUN: BRCAaccuTest PLUS, PANEL: BRCAaccuTest PLUS, DISEASE: HBOC, INSTRUMENT: Illumina MiSeqDx. The main content area has tabs for “VARIANTS”, “OVERVIEW”, and “REPORT”, with “VARIANTS” currently selected (yellow circle with 4). A sub-table for “SNV/Indel: BRCA1 missense_variant NM_007294.3 c.4837A>G p.S1613G” is shown, with “CNV” tab also visible. The table lists variants with columns for Position, Gene, Type, Consequence, NT Change, NT Change (BIC), and AA Change. A filter sidebar on the left shows “17/17(R:0)”, “Filter selection: All”, “Analysis actions: P, LP, US, LB, B, Report”, “Pathogenicity”, “Report”, “IGV”, and “False”. A download section at the bottom left says “Download a file” with icons for PDF, CSV, and VCF. The bottom right corner shows “System version: 1.4.10”.

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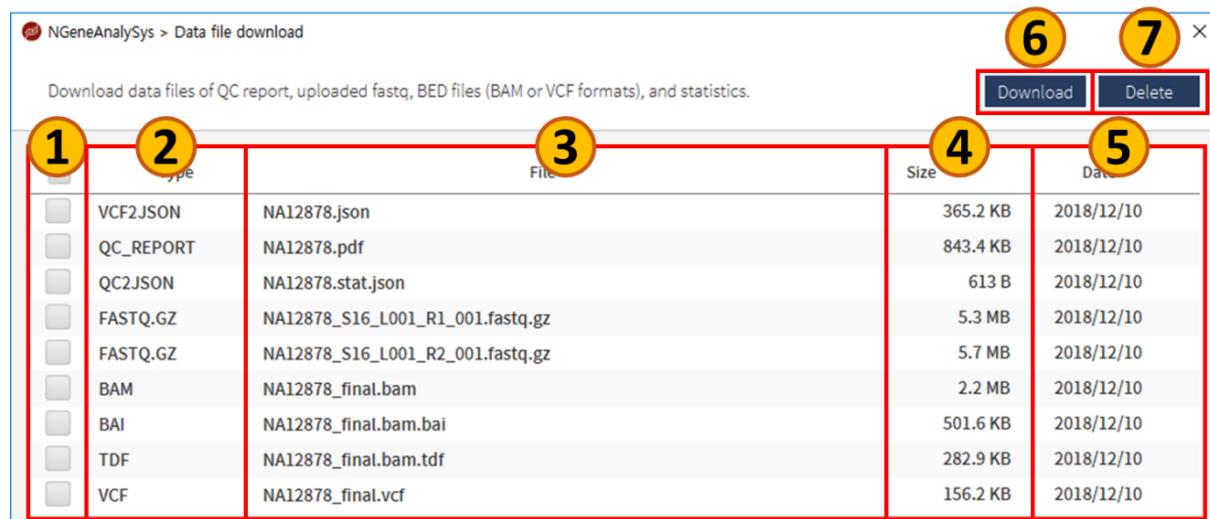Download data files of QC report, uploaded fastq, BED files (BAM or VCF formats), and statistics.				
Type	File	Size	Date	
<input type="checkbox"/> VCF2JSON	NA12878.json	365.2 KB	2018/12/10	
<input type="checkbox"/> QC_REPORT	NA12878.pdf	843.4 KB	2018/12/10	
<input type="checkbox"/> QC2JSON	NA12878.stat.json	613 B	2018/12/10	
<input type="checkbox"/> FASTQ.GZ	NA12878_S16_L001_R1_001.fastq.gz	5.3 MB	2018/12/10	
<input type="checkbox"/> FASTQ.GZ	NA12878_S16_L001_R2_001.fastq.gz	5.7 MB	2018/12/10	
<input type="checkbox"/> BAM	NA12878_final.bam	2.2 MB	2018/12/10	
<input type="checkbox"/> BAI	NA12878_final.bam.bai	501.6 KB	2018/12/10	
<input type="checkbox"/> TDF	NA12878_final.bam.tdf	282.9 KB	2018/12/10	
<input type="checkbox"/> VCF	NA12878_final.vcf	156.2 KB	2018/12/10	

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NA12878 RUN: BRCAaccuTest PLUS PANEL: BRCAaccuTest PLUS DISEASE: HBOC INSTRUMENT: Illumina MiSeqDx

Pa...	W...	Re...	Gene	Transcript...	Type	Consequence	NT Change	NT Change (BIC)	AA Ch...
B			BRCA1	NM_0072...	SNV	missense_variant	c.4837A>G	c.4956A>G	p.S16...
B			BRCA1	NM_0072...	Del	intron_variant	c.548-58delT	c.667-177delT	
B			BRCA1	NM_0072...	SNV	synonymous_variant	c.2082C>T	c.2201C>T	p.S69...
B			BRCA1	NM_0072...	SNV	synonymous_variant	c.2311T>C	c.2430T>C	p.L77...
B			BRCA1	NM_0072...	SNV	missense_variant	c.2612C>T	c.2731C>T	p.P87...
B			BRCA1	NM_0072...	SNV	missense_variant	c.3113A>G	c.3232A>G	p.E10...
B			BRCA1	NM_0072...	SNV	missense_variant	c.3548A>G	c.3667A>G	p.K11...
B			BRCA1	NM_0072...	SNV	synonymous_variant	c.4308T>C	c.4427T>C	p.S14...
B			BRCA2	NM_0000...	SNV	intron_variant	c.4485-63C>G	c.4604-182C>G	
B			BRCA2	NM_0000...	SNV	missense_variant	c.7397T>C	c.7625T>C	p.V24...
B			BRCA2	NM_0000...	SNV	synonymous_variant	c.7242A>G	c.7470A>G	p.S24...

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분류	아이콘	설명
Pathogenic	P	Causative for clinical phenotype.
Likely Pathogenic	LP	Probably causative for clinical phenotype.
Uncertain Significance	US	Unknown effect on clinical phenotype.
Likely Benign	LB	Probably not causative for a clinical phenotype.
Benign	B	Not causative for a clinical phenotype.

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The screenshot shows a software interface for variant analysis. On the left, there's a sidebar with filters for 'Pathogenicity' (P, LP, US, LB, B), 'Report', 'IGV', and 'False'. The main area displays a table of variants for the BRCA1 gene. The columns include: Position (17/17(R:1)), Filter selection (All), Analysis actions (P, LP, US, LB, B), Gene (BRCA1), Transcript ID (NM_007294.3), Type (highlighted with a red box), Consequence, NT Change, NT Change (BIC), and AA Change. A red box highlights the 'Type' column header. The table lists various variants, including some with 'R' (Risk) status and 'B' (Benign) pathogenicity.

SNV/Indel : BRCA1 intron_variant NM_007294.3 c.548-58delT		CNV		Reset Columns				
17/17(R:1)	Filter selection	All	New filter	Type	Consequence	NT Change	NT Change (BIC)	AA Change
R	BRCA1	NM_007294.3	c.548-58delT	Del	intron_variant	c.548-58delT	c.667-177delT	
B	BRCA1	NM_007294.3	c.4837A>G	SNV	missense_variant	c.4837A>G	c.4956A>G	p.S167L
B	BRCA1	NM_007294.3	c.2082C>T	SNV	synonymous_variant	c.2082C>T	c.2201C>T	p.S694S
B	BRCA1	NM_007294.3	c.2311T>C	SNV	synonymous_variant	c.2311T>C	c.2430T>C	p.L771L
B	BRCA1	NM_007294.3	c.2612C>T	SNV	missense_variant	c.2612C>T	c.2731C>T	p.P87H
B	BRCA1	NM_007294.3	c.3113A>G	SNV	missense_variant	c.3113A>G	c.3232A>G	p.E104Q
B	BRCA1	NM_007294.3	c.3548A>G	SNV	missense_variant	c.3548A>G	c.3667A>G	p.K119E
B	BRCA1	NM_007294.3	c.4308T>C	SNV	synonymous_variant	c.4308T>C	c.4427T>C	p.S143S
B	BRCA1	NM_007294.3	c.4485-63C>G	SNV	intron_variant	c.4485-63C>G	c.4604-182C>G	
B	BRCA2	NM_000003.3	c.7806-14T>C	SNV	intron_variant	c.7806-14T>C	c.8034-242T>C	
B	BRCA2	NM_000003.3	c.7397T>C	SNV	missense_variant	c.7397T>C	c.7625T>C	p.V249F

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

SNV/Indel Variant Filter 설정 화면 (Population Frequency)

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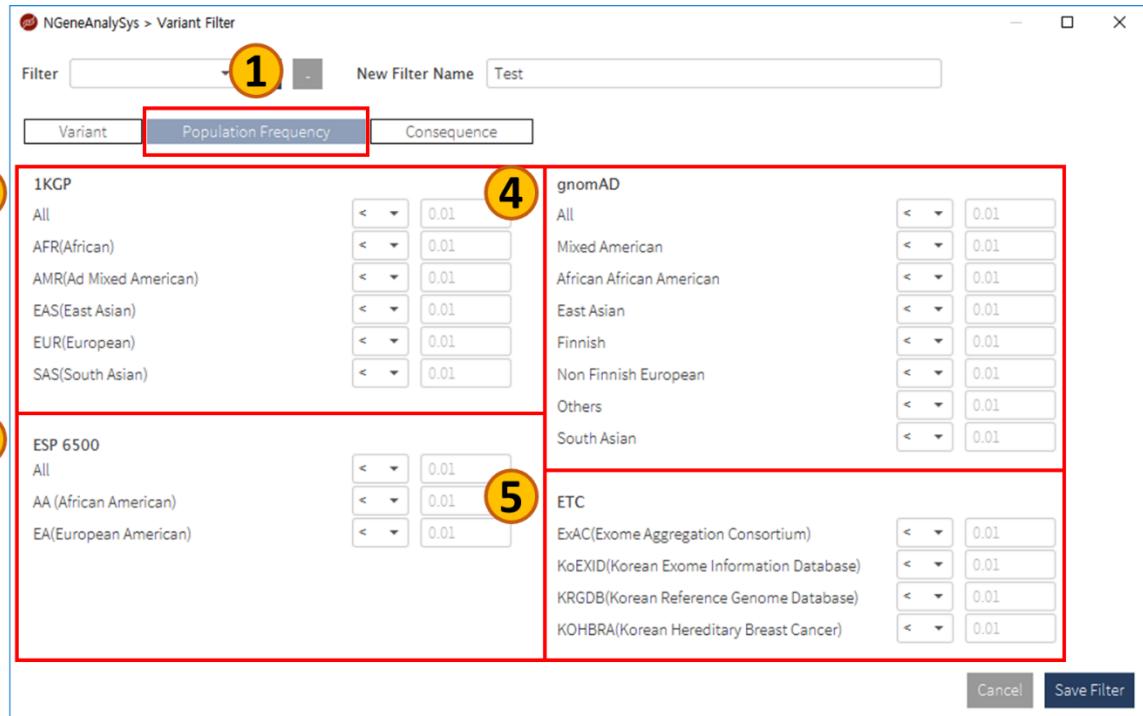SNV/Indel Variant Filter 설정 화면 (Consequence)

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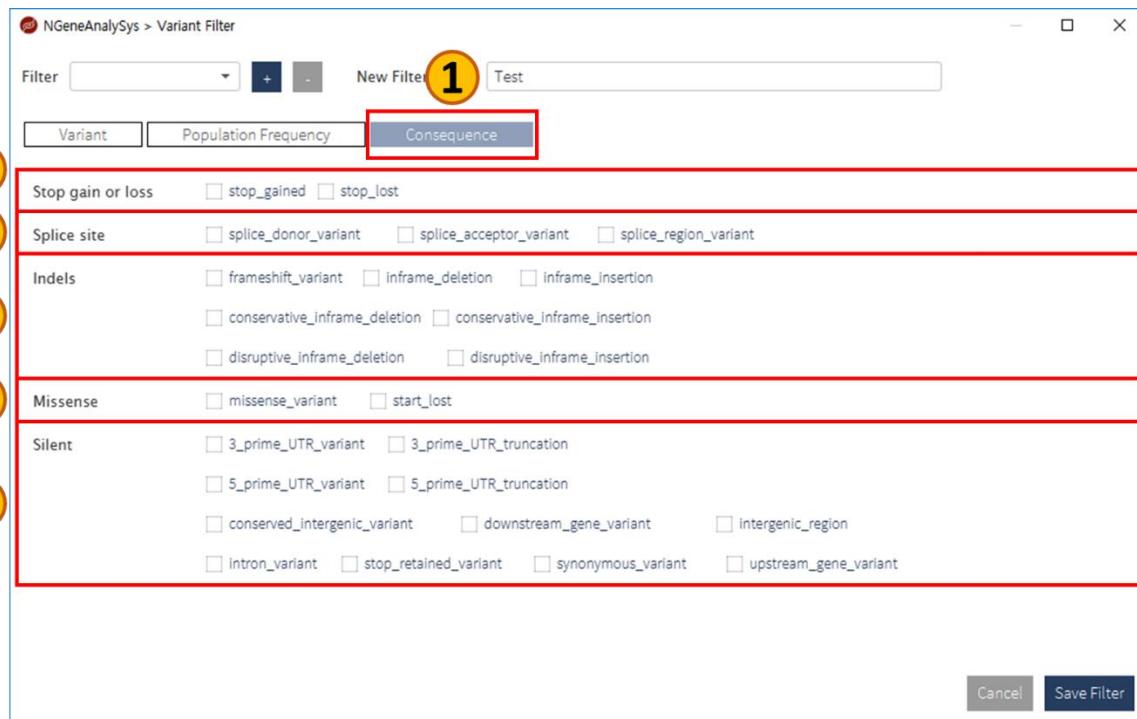
The screenshot shows the 'Variant Filter' configuration interface. At the top, there is a 'Filter' dropdown menu, a '+' button (highlighted with a red box), a '-' button, and a 'New Filter Name' input field containing 'Test' (also highlighted with a red box). Below these are several filter categories with checkboxes:

- Pathogenicity:** P(Pathogenic), LP(Likely Pathogenic), US(Uncertain Significance), LB(Likely Benign), B(Benign)
- ClinVar:** P(Pathogenic), LP(Likely Pathogenic), US(Uncertain Significance), LB(Likely Benign), B(Benign), Drug Response
- Variant Type:** SNV, Ins(Insertion), Del(Deletion), MNP(Multiple-nucleotide polymorphism), Complex(Multiple-nucleotides and an InDel)
- Chromosome:** chr7, chr9
- Gene:** FLT3, CEBPA
- dbSNP:** dbSNP ID
- COSMIC:** COSMIC ID, COSMIC Occurrence (Tissue selection dropdown), Allele Fraction (30 ~ 50%), Depth (30 ~ 50%), Alt Count (30 ~ 50%)
- Reads Coverage:** Allele Fraction (30 ~ 50%), Depth (30 ~ 50%), Alt Count (30 ~ 50%)
- Warning(Low Confidence):** All, Default, Uncheck

At the bottom right of the window are two buttons: 'Cancel' and 'Save Filter' (highlighted with a red box).

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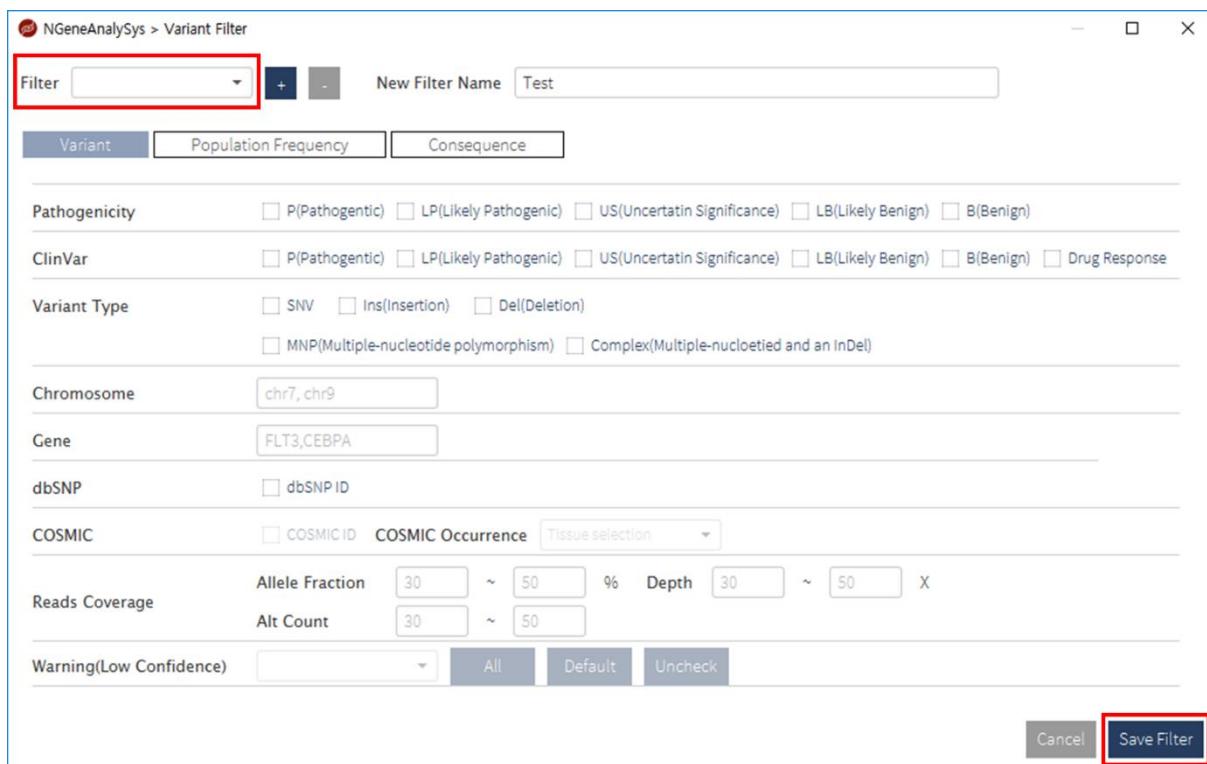
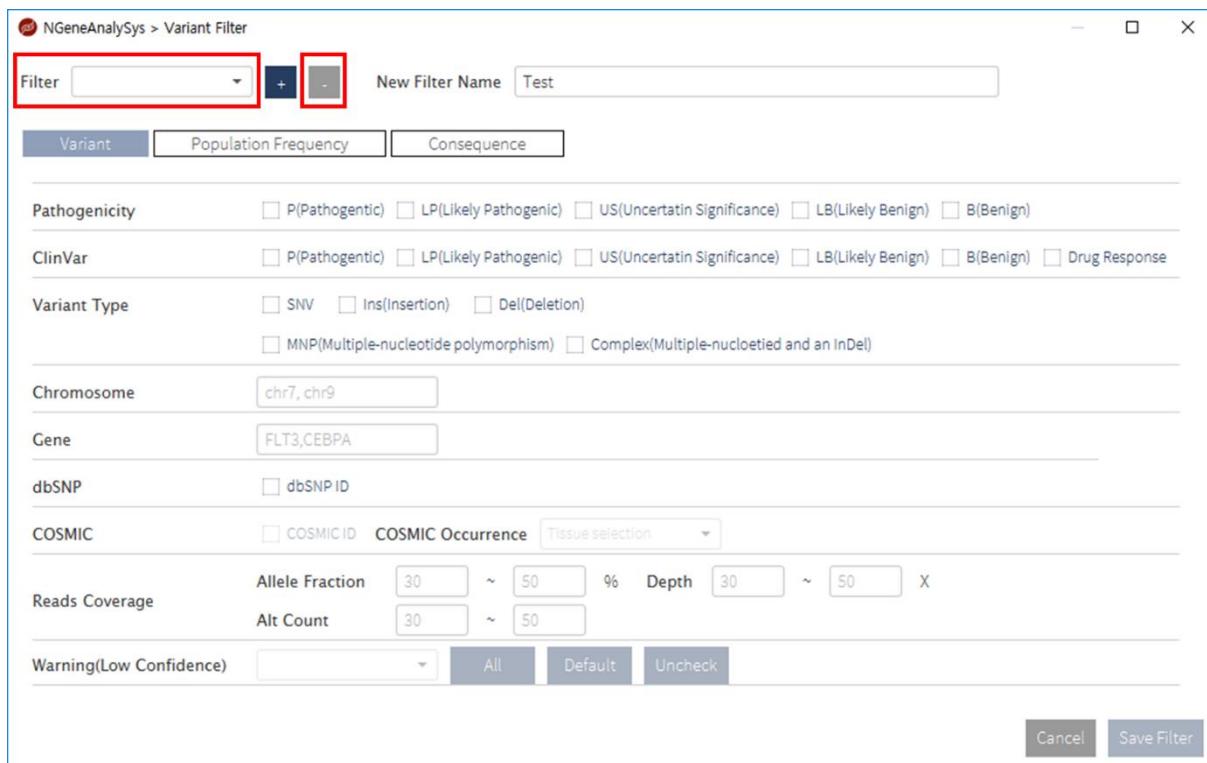
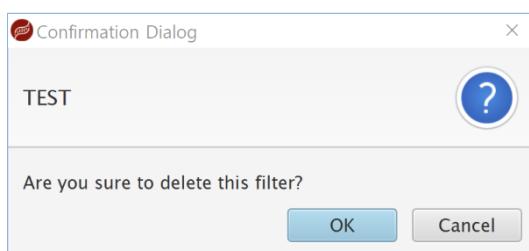
The screenshot shows the NGeneAnalySys software interface. At the top, there are tabs for DASHBOARD, RESULTS, and SAMPLES (which is currently selected). Below the tabs, it displays test details: RUN: BRCAaccuTest PLUS, PANEL: BRCAaccuTest PLUS v1, DISEASE: HBOC, SEQUENCER: Illumina MiSeqDx. In the center, under the 'SAMPLES' tab, there's a list of variants. One specific variant is highlighted: SNV/Indel : BRCA1 frameshift_variant NM_007300.3 c.5329dupC p.Q1777fs. A red box surrounds the 'Variant Detail' button in the dropdown menu that appears when this item is selected. The bottom right corner of the window shows 'System version: 14.0.0'.

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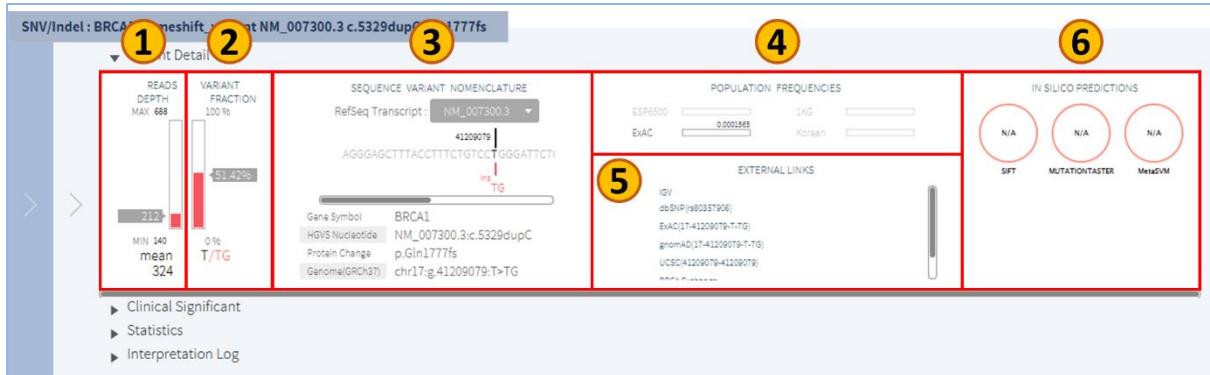The screenshot shows the 'Statistics' section of the SNV/Indel variant details. It includes three circular icons labeled 1, 2, and 3, each representing a different context for variant interpretation:

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The screenshot shows the 'Interpretation Log' section of the SNV/Indel variant details. It displays a table with the following columns:

Date	Type	User	Previous	Current	Comment
No content in table					

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
Pathogenicity 일괄 변경

The screenshot shows a software interface for variant analysis. On the left, a sidebar titled 'SNV/Indel: BRCA1 intron_variant NM_007294.3 c.548-58delT' displays '17/17(R:1)' variants. It includes a 'Filter selection' dropdown set to 'All', a 'New filter' button, and an 'Analysis actions' section with checkboxes for P, LP, US, LB, and B, with 'B' selected. A red box highlights the 'Pathogenicity' button. Below these are 'Report', 'IGV', and 'False' buttons, and download options for 'Download a file' (FASTA, VCF, CSV). On the right, a main table titled 'CNV' lists variants. The columns are: Pa..., W..., Re..., Gene, Transcript..., Type, Consequence, NT Change, NT Change (BIC), AA Ch. The first row shows a variant with 'B' status, 'BRCA1' gene, 'NM_0072...' transcript, 'Del' type, 'intron_variant' consequence, 'c.548-58delT' NT change, 'c.667-177delT' BIC, and 'p.S16...' AA change.

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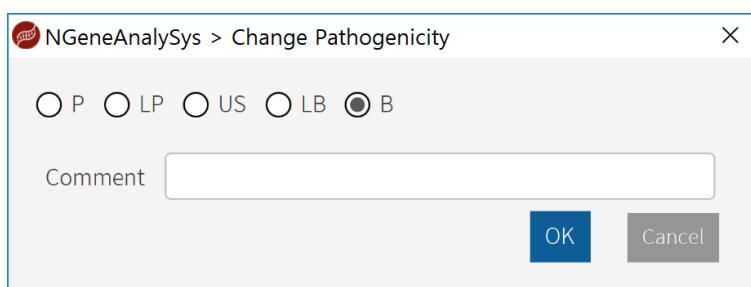
The screenshot shows the software's main window for a specific variant. Section 1 (Variant Details) includes 'Variant Details' and 'Clinical Significance'. Section 2 (Predicted Pathogenicity) contains 'Predicted Pathogenicity' and 'Change Pathogenicity' buttons, both highlighted with a red box. Section 3 (Pathogenicity classification) lists 'PREDICTION' (BIC, CUNIVAR, ENIGMA), 'CLINVAR', and 'ENIGMA' with their respective pathogenicity scores. Section 4 (Detailed Comment) shows a large box for 'Add to Report' containing the 'PVS1' classification with its definition: 'Null variant (nonsense, frameshift, canonical +1 or 2 splice sites, initiation codon, single or multixon deletion) in a gene where LOF is a known mechanism of disease. (Loss of function is a common mechanism of action in this gene (BRCA1). This variant'.

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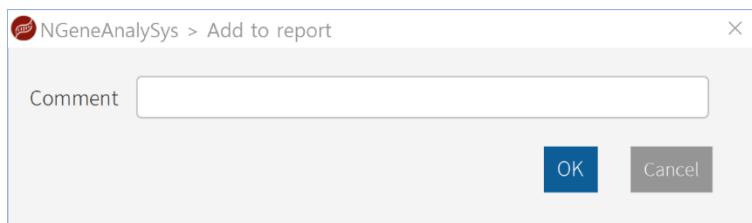
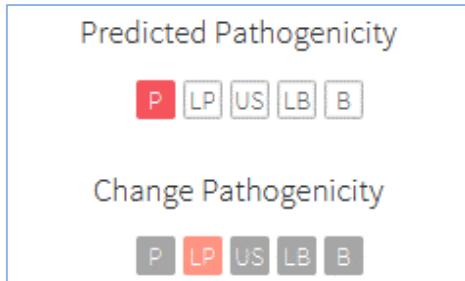
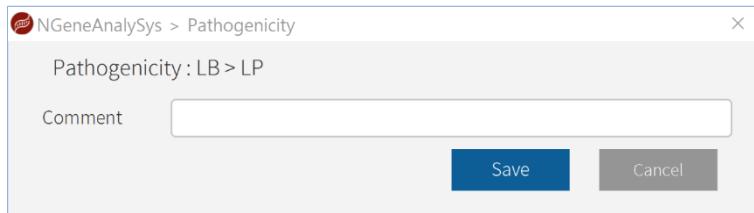
The screenshot shows a variant analysis interface with a sidebar titled 'Analysis actions'. In the sidebar, the 'Tier' button is highlighted with a red box. The main area displays a table of variants with columns for Tier, Gene, Type, Consequence, NT Change, AA Change (Single), Chr, and Start Pos.

	Tier	Gene	Type	Consequence	NT Change	AA Change (Single)	Chr	Start Pos
T1	FLT3	Ins	conservative_infram...	c.1786_1806dupGAAT...	p.K602_W603insEY...	chr13	2860	
T1	FLT3	Ins	disruptive_infram...	c.1805_1806insGCGCT...	p.K602_W603insR...	chr13	2860	
T1	FLT3	Ins	disruptive_infram...	c.1719_1724dupAAGCCA	p.Q575_L576insSQ	chr13	2860	
T1	FLT3	SNV	missense_variant	c.2503G>T	p.D835Y	chr13	2853	
T1	RUNX1	Del	conservative_infram...	c.331_336delACCTTG	p.T111_L112del	chr21	3623	
T1	RUNX1	Ins	frameshift_variant	c.1136_1137insA	p.Y380fs	chr21	3610	
T2	DNMT3A	SNV	missense_variant	c.2645G>A	p.R882H	chr2	2543	
T3	ANKRD26	Del	splice_region_vari...	c.1986_4delT		chr10	2733	
T3	ASXL1	Del	conservative_infram...	c.25_27delAAG	p.K9del	chr20	3094	
T3	BCR	SNV	missense_variant	c.331C>G	p.A1104G	chr22	2363	
T3	CBL	Del	disruptive_infram...	c.125_127delACC	p.H42del	chr11	1190	

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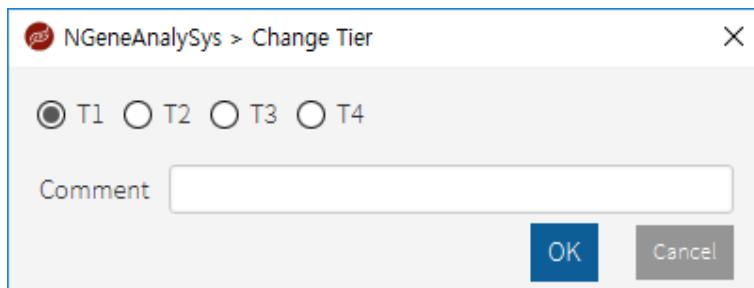Variant Detail
Interpretation

Tier 1 2 3 4 5 6 7

Type	Evidence	Evidence Comment	Primary	Delete
Prognosis	C	Poor in MDS	<input type="radio"/>	<input type="checkbox"/>
Prognosis	A	Adverse in AML (ELN 2017), Poor in AML (NCCN 2017)	<input type="radio"/>	<input type="checkbox"/>
Diagnosis	A	Suggestive AML, MDS (NCCN 2018)	<input checked="" type="radio"/>	<input type="checkbox"/>
Therapeutic	D	ROS045337 (phase 1, NCT01677780)	<input type="radio"/>	<input type="checkbox"/>

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	Gene	Exon	CNV	Report	Warning	Domain	Total Amplicon	Number of Amplicon Copy			Comment
								Loss	Normal	Gain	
	BRCA1	2	Copy Loss		R	RING	1	1	0	0	2019-03-07, Report(N->Y), Test
	BRCA1	3	Copy Loss		R	RING	1	1	0	0	2019-03-07, Report(N->Y), Test
	BRCA1	4	Copy Loss		R	RING	1	1	0	0	2019-03-07, Report(N->Y), Test
	BRCA1	5	Copy Loss		R		1	1	0	0	
	BRCA1	6	Copy Loss		R		1	1	0	0	
	BRCA1	7	Copy Loss		R		1	1	0	0	
	BRCA1	8	Copy Loss		R		1	1	0	0	
	BRCA1	9	Copy Loss		R		1	1	0	0	
	BRCA1	10	Copy Loss		R		26	23	3	0	
	BRCA1	11	Copy Loss		R	PALB2	1	1	0	0	
	BRCA1	12	Copy Loss		R	PALB2	2	2	0	0	
	BRCA1	13	Copy Loss		R		2	2	0	0	
	BRCA1	14	Copy Loss		R		2	2	0	0	
	BRCA1	15	Copy Loss		R	BRCT1	3	3	0	0	
	BRCA1	16	Copy Loss		R	BRCT1	1	1	0	0	

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Clinical Report	
Test Information	Client Information
Test No.	Ordering Organization
Test Name	Ordering Physician
Contact:	
Specimens	Patient
Specimen Type: 血液	Name: 患者名 Date of Birth: 生年月日 Received Date: 受取日 Patient ID: 患者ID Gender: 性別
Summary	
Positive clinically significant mutation identified: BRCA2, NM_000093.2(c.146C>T, p.(Q49X)), hemizygote	
Result/Conclusion	
Details, BRCA2	
N/C Change (SNP/Frame): N/C Change ¹ (SNP/Frame): AA Change ² Zygosity: Type: Classification	N/C Change (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) N/C Change ¹ (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) AA Change ² : NM_000093.2(c.146C>T, p.(Q49X)) Zygosity: hemizygote Type: SNP B Classification: L
N/C Change (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) N/C Change ¹ (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) AA Change ² : NM_000093.2(c.146C>T, p.(Q49X)) Zygosity: hemizygote Type: SNP B Classification: L	N/C Change (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) N/C Change ¹ (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) AA Change ² : NM_000093.2(c.146C>T, p.(Q49X)) Zygosity: hemizygote Type: SNP B Classification: L
Test Information	
1. Specimen: 血液 2. Clinical information 3. Test method	
Details, BRCA2	
N/C Change (SNP/Frame): N/C Change ¹ (SNP/Frame): AA Change ² Zygosity: Type: Classification	N/C Change (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) N/C Change ¹ (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) AA Change ² : NM_000093.2(c.146C>T, p.(Q49X)) Zygosity: hemizygote Type: SNP B Classification: L
N/C Change (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) N/C Change ¹ (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) AA Change ² : NM_000093.2(c.146C>T, p.(Q49X)) Zygosity: hemizygote Type: SNP B Classification: L	N/C Change (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) N/C Change ¹ (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) AA Change ² : NM_000093.2(c.146C>T, p.(Q49X)) Zygosity: hemizygote Type: SNP B Classification: L
Test Information	
1. Specimen: 血液 2. Clinical information 3. Test method	
Details, BRCA2	
N/C Change (SNP/Frame): N/C Change ¹ (SNP/Frame): AA Change ² Zygosity: Type: Classification	N/C Change (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) N/C Change ¹ (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) AA Change ² : NM_000093.2(c.146C>T, p.(Q49X)) Zygosity: hemizygote Type: SNP B Classification: L
N/C Change (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) N/C Change ¹ (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) AA Change ² : NM_000093.2(c.146C>T, p.(Q49X)) Zygosity: hemizygote Type: SNP B Classification: L	N/C Change (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) N/C Change ¹ (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) AA Change ² : NM_000093.2(c.146C>T, p.(Q49X)) Zygosity: hemizygote Type: SNP B Classification: L
Test Information	
1. Specimen: 血液 2. Clinical information 3. Test method	

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
Clinical Report	
2018-05-16	
Test Information	Client Information
Test No.	Ordering Organization
Test Name	Ordering Physician
Contact:	
Specimens	Patient
Specimen Type: 血液	Name: 患者名 Date of Birth: 生年月日 Received Date: 受取日 Patient ID: 患者ID Gender: 性別
Summary	
Positive clinically significant mutation identified: BRCA2, NM_000093.2(c.146C>T, p.(Q49X)), hemizygote	
Result/Conclusion	
Details, BRCA2	
N/C Change (SNP/Frame): N/C Change ¹ (SNP/Frame): AA Change ² Zygosity: Type: Classification	N/C Change (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) N/C Change ¹ (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) AA Change ² : NM_000093.2(c.146C>T, p.(Q49X)) Zygosity: hemizygote Type: SNP B Classification: L
N/C Change (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) N/C Change ¹ (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) AA Change ² : NM_000093.2(c.146C>T, p.(Q49X)) Zygosity: hemizygote Type: SNP B Classification: L	N/C Change (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) N/C Change ¹ (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) AA Change ² : NM_000093.2(c.146C>T, p.(Q49X)) Zygosity: hemizygote Type: SNP B Classification: L
Test Information	
1. Specimen: 血液	
Details, BRCA2	
N/C Change (SNP/Frame): N/C Change ¹ (SNP/Frame): AA Change ² Zygosity: Type: Classification	N/C Change (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) N/C Change ¹ (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) AA Change ² : NM_000093.2(c.146C>T, p.(Q49X)) Zygosity: hemizygote Type: SNP B Classification: L
N/C Change (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) N/C Change ¹ (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) AA Change ² : NM_000093.2(c.146C>T, p.(Q49X)) Zygosity: hemizygote Type: SNP B Classification: L	N/C Change (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) N/C Change ¹ (SNP/Frame): NM_000093.2(c.146C>T, p.(Q49X)) AA Change ² : NM_000093.2(c.146C>T, p.(Q49X)) Zygosity: hemizygote Type: SNP B Classification: L
Test Information	
1. Specimen: 血液	

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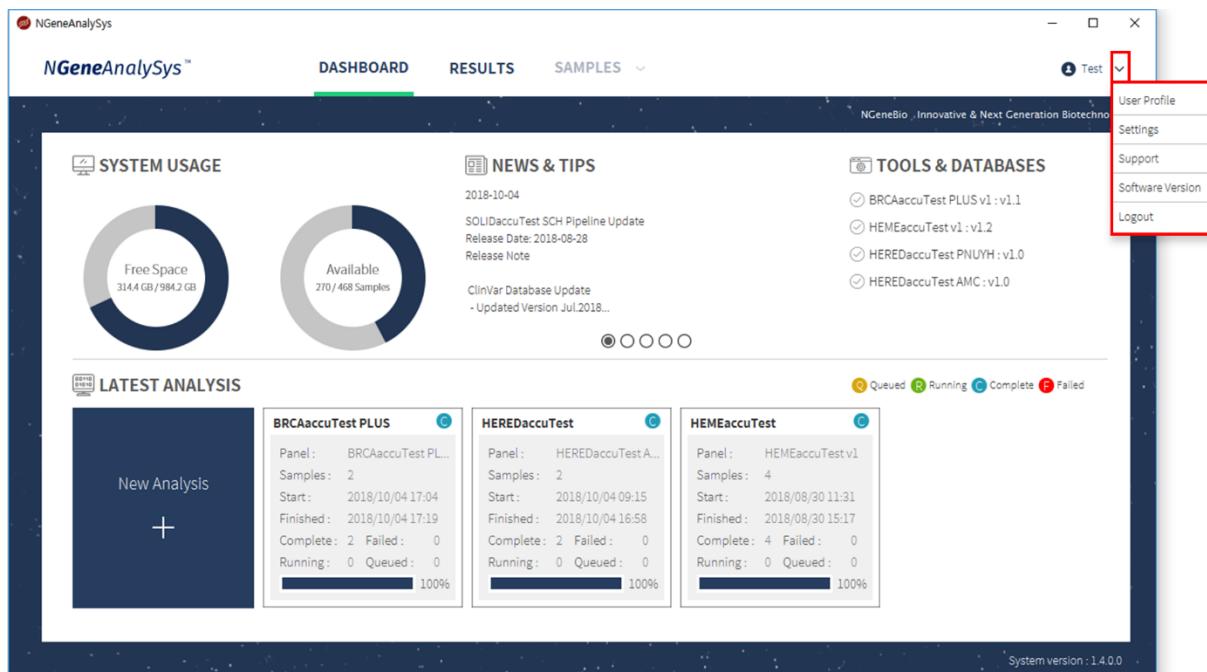The screenshot shows the 'User Profile' window with two main sections highlighted by a red box: 'Change Password' and 'User Information'. The 'Change Password' section contains fields for Current Password, New Password, and New Confirm Password, with a note below about password requirements. The 'User Information' section contains fields for Organization, Department, Address, Phone, and E-mail. At the bottom right are 'Save' and 'Close' buttons, with a yellow circle labeled '3' over the 'Save' button.

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The screenshot shows the NGeneAnalySys software interface. At the top, there is a navigation bar with tabs: DASHBOARD, RESULTS, SAMPLES, and MANAGER (which is highlighted with a red box and a yellow circle containing the number 1). Below the navigation bar, there is a sub-menu bar with several items: Analysis Status, User Account, Panel, Report, News&Tips, System Logs, and Clinical Variant DB. A yellow circle with the number 2 is placed over the User Account item. The main content area displays a table of analysis runs. The table has columns: Run ID, Run Name, Status, Login ID, User, Group, Date, Restart, and Delete. There are three rows of data:

Run ID	Run Name	Status	Login ID	User	Group	Date	Restart	Delete
3	BRCA_PLUS_FFPE	COMPLETE	admin	Administrator	Administrators	2018-05-14 14:34:23		
2	BRCA_PLUS_BLOOD	COMPLETE	admin	Administrator	Administrators	2018-05-14 14:34:07		
1	BRCA	COMPLETE	admin	Administrator	Administrators	2018-05-14 14:33:54		

A yellow circle with the number 3 is placed over the first row of the table. At the bottom of the table, there are navigation buttons for page selection (Page, GO) and a Refresh button. The bottom right corner of the interface shows the text "System version : 1.4.8.0".

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The screenshot shows the 'Analysis Status' section of the page. A red box highlights the 'Analysis Status' button at the top left. A yellow circle labeled '1' is over this button. Another yellow circle labeled '2' is over the table rows. A third yellow circle labeled '3' is over the 'Refresh' button in the top right corner of the table area.

Run ID	Run Name	Status	Login ID	User	Group	Date	Restart	Delete
3	BRCA_PLUS_FFPE	COMPLETE	admin	Administrator	Administrators	2018-05-14 14:34:23		
2	BRCA_PLUS_BLOOD	COMPLETE	admin	Administrator	Administrators	2018-05-14 14:34:07		
1	BRCA	COMPLETE	admin	Administrator	Administrators	2018-05-14 14:33:54		

Below the table are navigation controls: '<<', '<', '1', '>', '>>', 'Page []', and 'GO'.

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The screenshot shows the 'User Account' tab selected in the top navigation bar. Below the tabs, there are two tables: 'Group' on the left and 'User' on the right. The 'User' table is highlighted with a red box and has three numbered callouts: 1 points to the 'User Account' tab; 2 points to the 'Group' table; and 3 points to the 'User' table.

ID	Name	Deleted	Modify
2	Test	N	
1	Administrators	N	

Login ID	Name	User Type	Group	Created At	Updated At	Deleted	Modify
Test	Test	RESEARC...	Test	2018-06-04 04:38..	-	N	
admin	Adminis...	ADMIN	Administrators	2018-06-04 00:44..	-	N	

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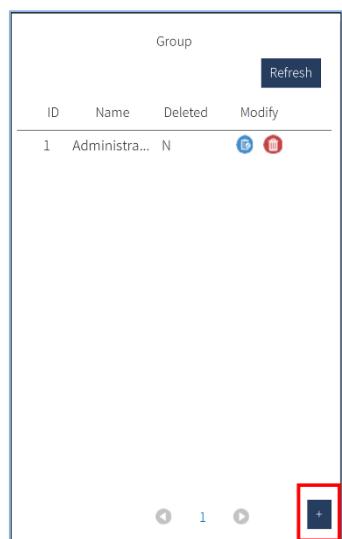
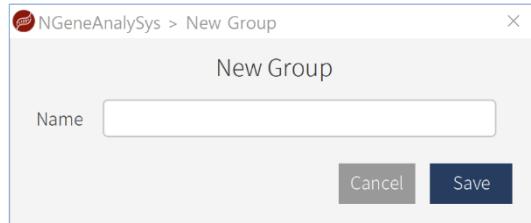
Group			
ID	Name	Deleted	Modify
2	Test	N	
1	Administrat...	N	

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
ID	Name	Deleted	Modify
2	Test	N	 
1	Administrator	N	 

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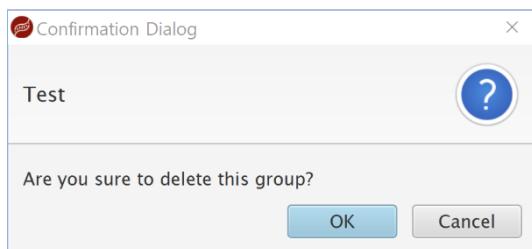The screenshot shows a table with one row of data. The columns are: Login ID, Name, User Type, Group, Created At, Updated At, Deleted, and Modify. The data in the first row is: admin, Adminis..., ADMIN, Administrators, 2018-05-17 08:26... - N. There is a 'Refresh' button at the top right and a 'Modify' button in the last column. A red box highlights the entire row, and a yellow circle with the number 2 is on its left. A yellow circle with the number 1 is on the 'Modify' button. A yellow circle with the number 3 is on the 'GO' button in the footer.

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The screenshot shows the same table structure as the previous one, but the 'Modify' button in the first row is now highlighted with a red box. A yellow circle with the number 1 is on the 'Modify' button. A yellow circle with the number 2 is on the 'GO' button in the footer. The 'Create' button in the bottom right corner is also highlighted with a red box, and a yellow circle with the number 3 is on it.

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The screenshot shows a 'New User' creation form. The fields filled in are:

- Login ID (*): Test
- Name (*): Test
- Password (*): [REDACTED]
- Confirm Password (*): [REDACTED]
- User Type (*): RESEARCHER
- User Group (*): Test
- Organization: [REDACTED]
- Department: [REDACTED]
- Address: [REDACTED]
- Phone: [REDACTED]
- E-mail: [REDACTED]

At the bottom right are two buttons: **Cancel** and **Save**, where **Save** is highlighted.

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User							
Login ID	Name	User Type	Group	Created At	Updated At	Deleted	Modify
Test	Test	RESEARC...	Test	2018-06-04 04:38...	-	N	 
admin	Adminis...	ADMIN	Administrators	2018-06-04 00:44...	-	N	 

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 NGeneAnalySys > Edit User

Edit User


Login ID (*)	Test
Name (*)	Test
Passsword (*)	
Confirm Passsword (*)	
User Type(*)	RESEARCHER
User Group(*)	Test
Organization	
Department	
Address	
Phone	
E-mail	

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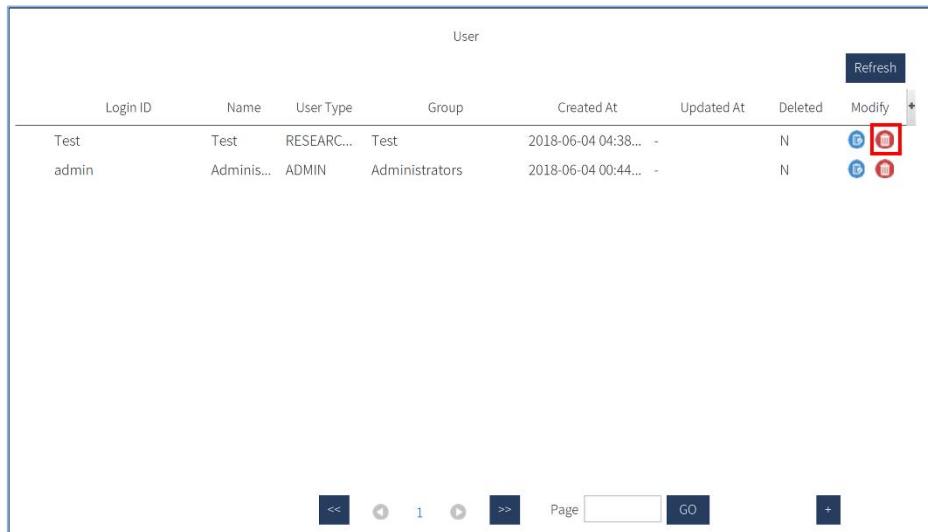
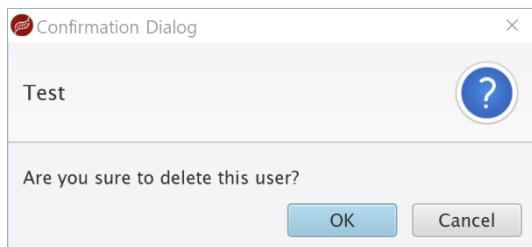
CNV for BRCAaccuTest Series Panel 화면

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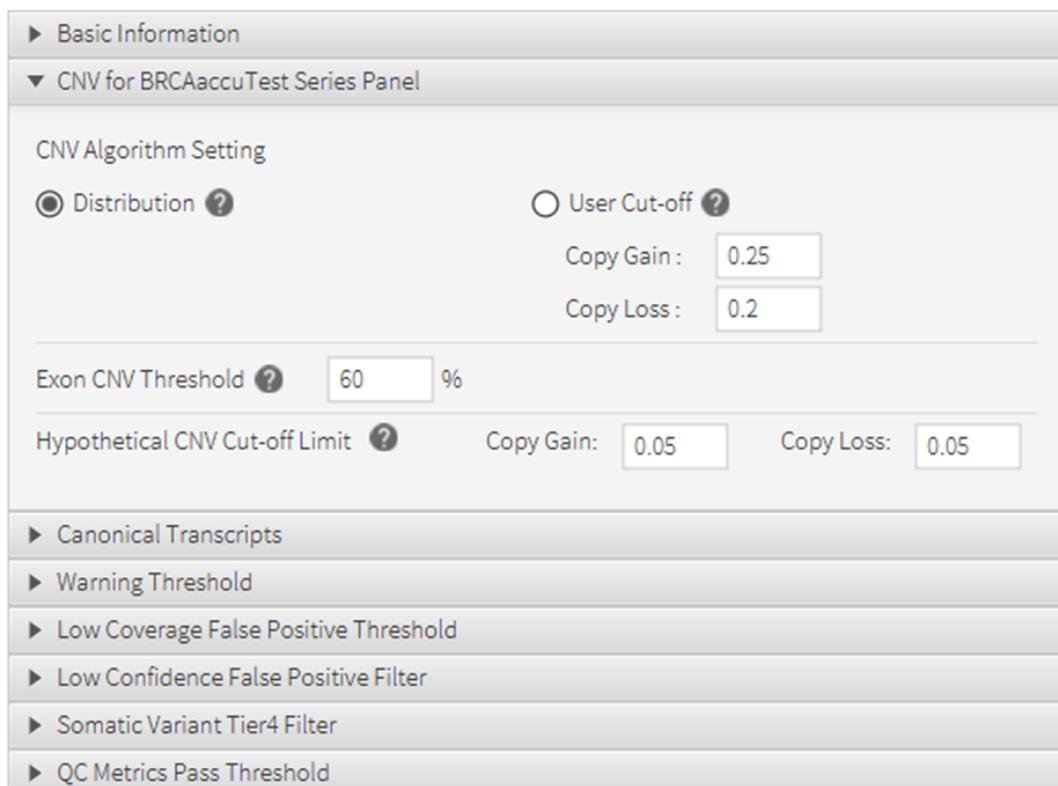
이름	설명
Distribution	CNV analysis using normal distribution
User Cut-off	CNV analysis applied with cut-off value(s) set by user (default Copy Gain = 0.25, Copy Loss = 0.2)
Exon CNV Threshold	Exon CNV determination according to the ratio of amplicon copy number (default = 60%)
Hypothetical CNV Cut-off Limit	Cut-off limit value for hypothetical gain or loss each amplicon (default = 0.05)

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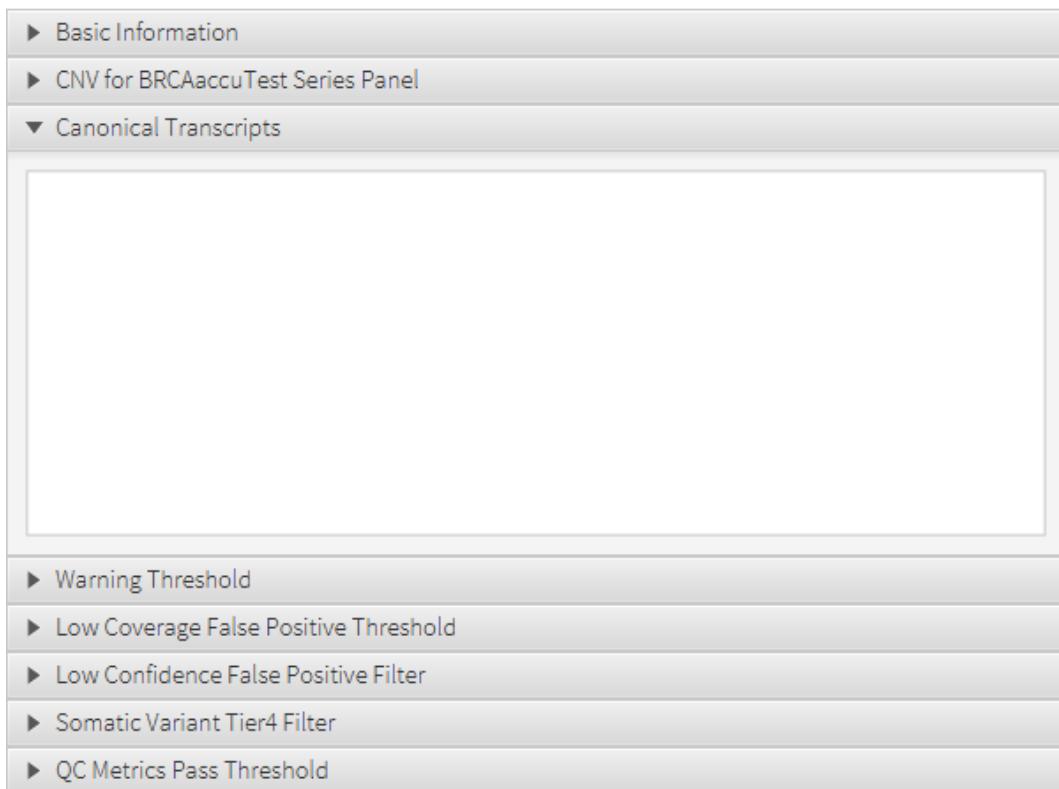
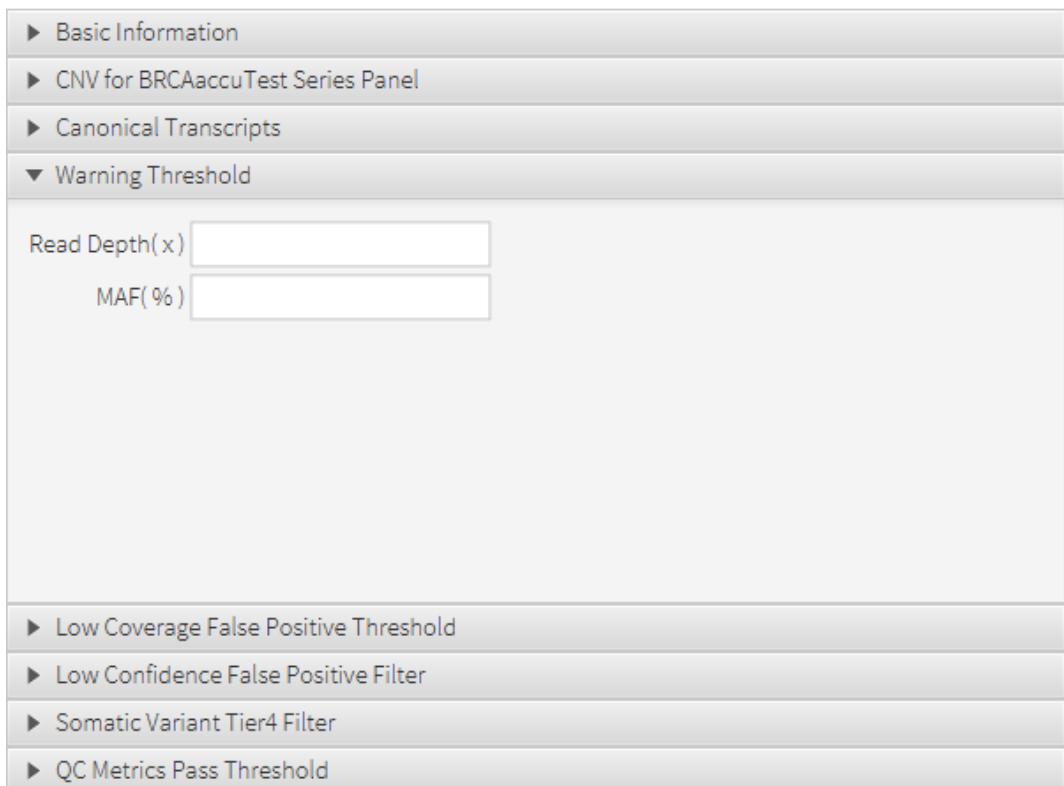
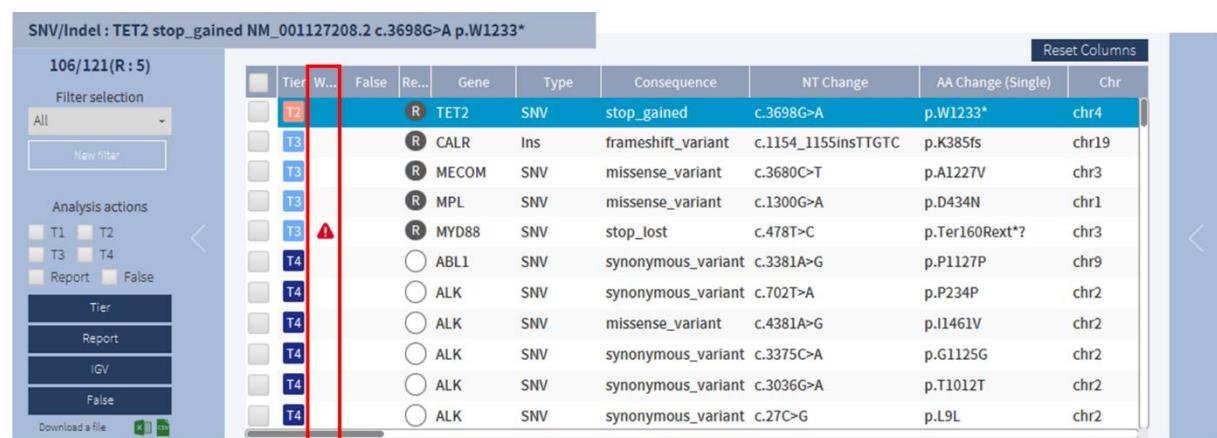
Low Coverage False Positive Threshold 화면

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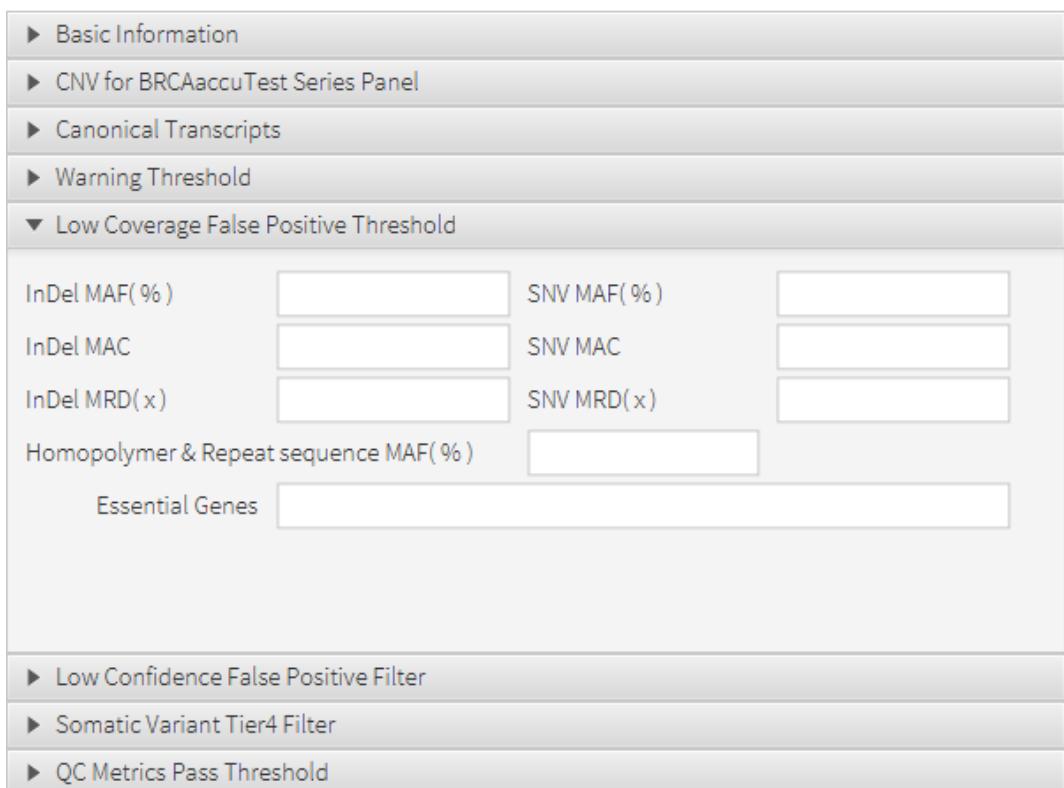▶ Canonical Transcripts

▶ Warning Threshold

▼ Low Coverage False Positive Threshold

InDel MAF(%)	<input type="text"/>	SNV MAF(%)	<input type="text"/>
InDel MAC	<input type="text"/>	SNV MAC	<input type="text"/>
InDel MRD(x)	<input type="text"/>	SNV MRD(x)	<input type="text"/>
Homopolymer & Repeat sequence MAF(%)		<input type="text"/>	
Essential Genes <input type="text"/>			

▶ Low Confidence False Positive Filter

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Low Confidence False Positive Filter 화면

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Low Confidence

- mapping_quality
- strand_artifact
- panel_of_normal
- fragment_length
- orientation_bias
- sequencing_error
- mapping_error
- snp_candidate
- t_lod

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▶ Canonical Transcripts

▶ Warning Threshold

▶ Low Coverage False Positive Threshold

▶ Low Confidence False Positive Filter

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
Below the table are three buttons: 'Save current setting to file', 'Select ROI file', and 'Save'.

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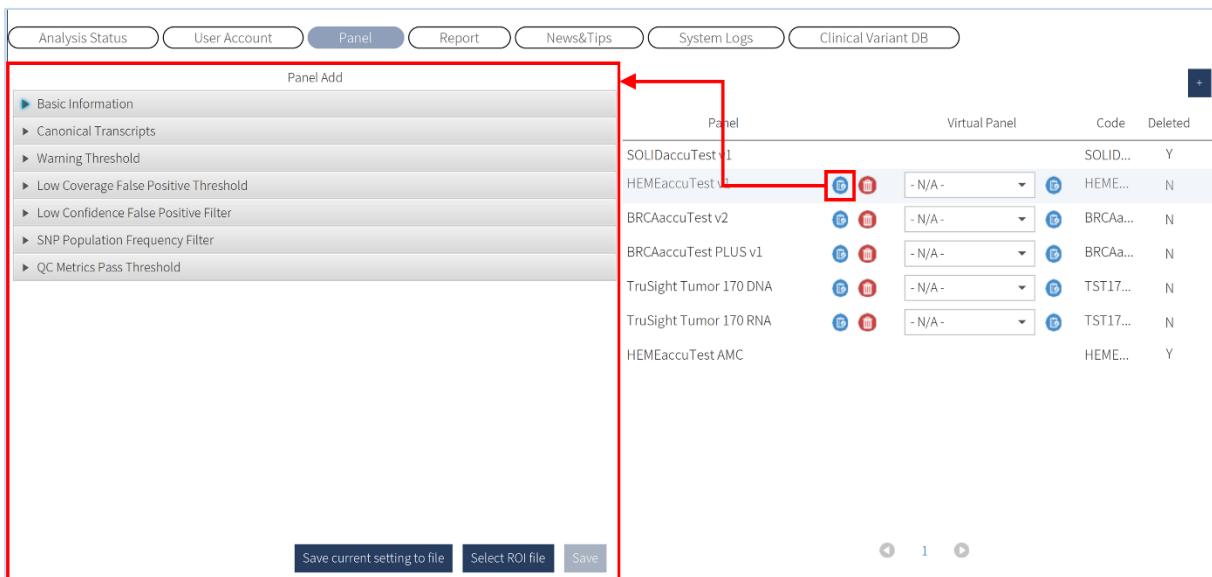
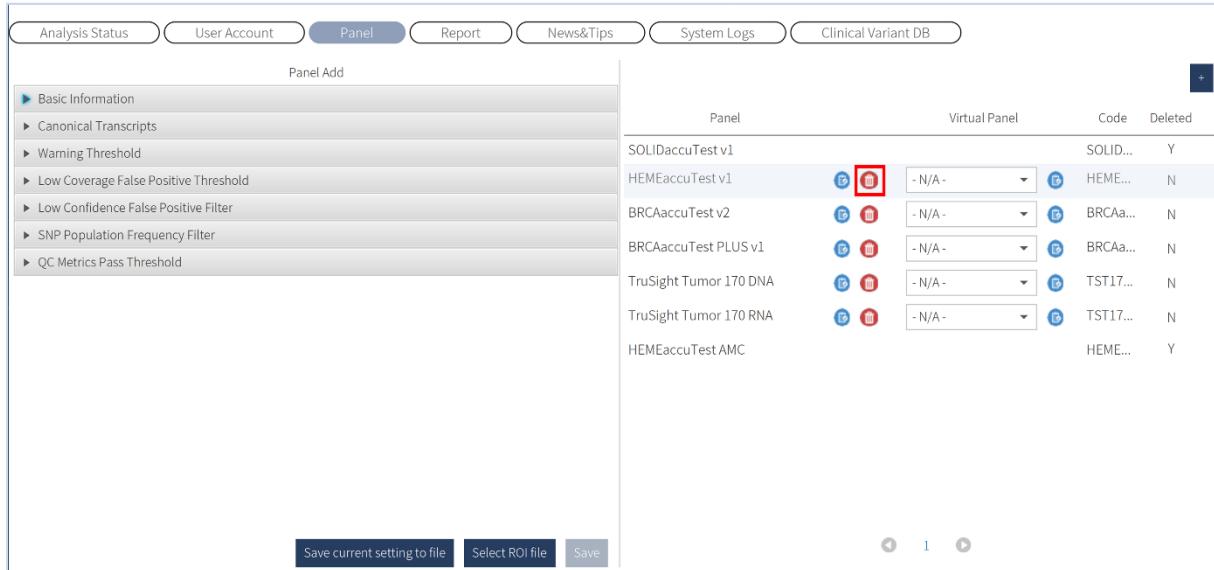
The screenshot shows a software interface for managing panels. At the top, there is a navigation bar with tabs: Analysis Status, User Account, Panel (which is selected and highlighted in blue), Report, News&Tips, System Logs, and Clinical Variant DB. Below the navigation bar is a sidebar titled "Panel Add" containing several expandable sections: Basic Information, Canonical Transcripts, Warning Threshold, Low Coverage False Positive Threshold, Low Confidence False Positive Filter, SNP Population Frequency Filter, and QC Metrics Pass Threshold. To the right of the sidebar is a main table area displaying a list of panels. The table has columns for "Panel", "Virtual Panel", "Code", and "Deleted". A "+" button is located in the top right corner of the table header. The data in the table is as follows:

Panel	Virtual Panel	Code	Deleted	
SOLIDaccuTest v1	- N/A -	SOLID...	Y	
HEMEaccuTest v1		- N/A -	HEME...	N
BRCAaccuTest v2		- N/A -	BRCAa...	N
BRCAaccuTest PLUS v1		- N/A -	BRCAa...	N
TruSight Tumor 170 DNA		- N/A -	TST17...	N
TruSight Tumor 170 RNA		- N/A -	TST17...	N
HEMEaccuTest AMC			HEME...	Y


At the bottom of the table area are three buttons: "Save current setting to file", "Select ROI file", and "Save". There are also navigation icons for back, forward, and search.

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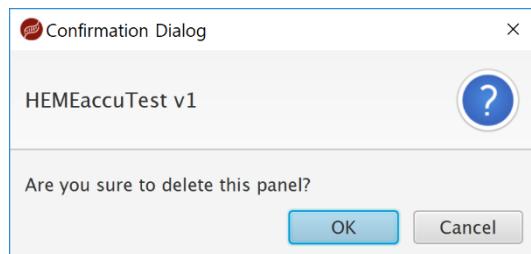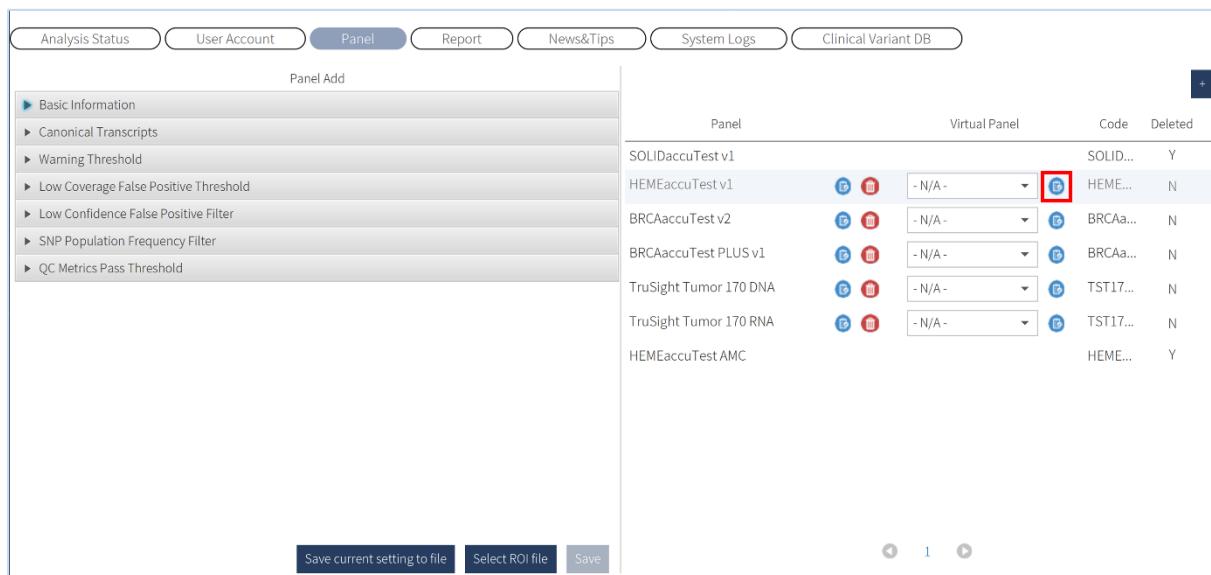The screenshot shows the 'Panel Add' section of a software interface. On the left, there's a sidebar with 'Basic Information' and several filter options like 'Canonical Transcripts', 'Warning Threshold', etc. The main area lists various panels with columns for 'Panel', 'Virtual Panel', 'Code', and 'Deleted'. A specific row for 'HEMEaccuTest v1' has its 'Virtual Panel' dropdown open, with the edit icon circled in red. At the bottom, there are buttons for 'Save current setting to file', 'Select ROI file', and 'Save'.

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
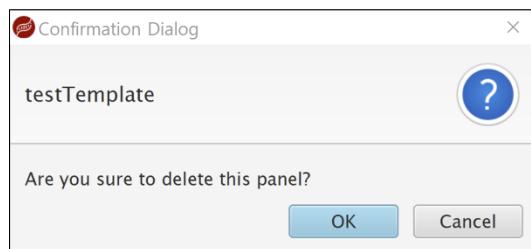
The screenshot shows a web-based application interface for managing news. At the top, there is a navigation bar with several tabs: Analysis Status, User Account, Panel, Report, News&Tips (which is highlighted with a red box and a yellow circle containing the number 1), System Logs, and Clinical Variant DB. Below the navigation bar, there is a list of news items. The first item in the list has its details displayed on the left side of the screen. This detailed view includes fields for Title (containing 'Test'), ID (containing '1'), and other metadata like Create Date (2018-05-17) and Update Date (2018-05-17). There are also 'Modify' and 'Delete' buttons. A red box highlights the 'News&Tips' tab and the detailed view area. Numbered circles (2, 3, 4) point to specific elements: 2 points to the 'Title' field in the detailed view; 3 points to the 'ID' field in the list; and 4 points to the 'Modify' button in the detailed view.

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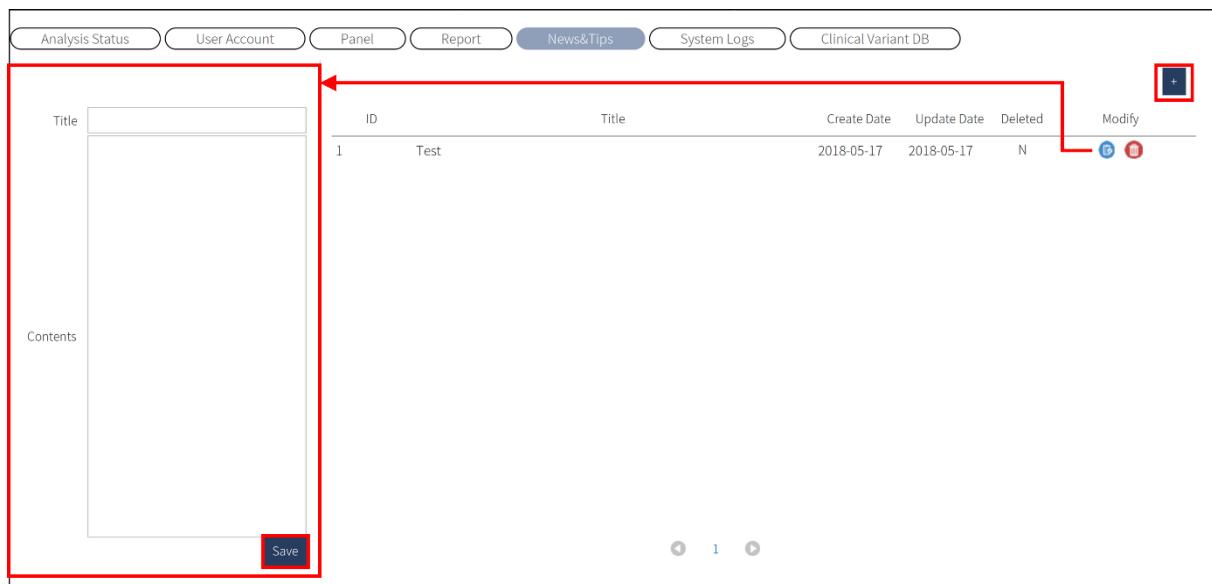
The screenshot shows a software interface for managing news and tips. At the top, there are several tabs: Analysis Status, User Account, Panel, Report, News&Tips (which is selected and highlighted in blue), System Logs, and Clinical Variant DB. Below the tabs, there's a search bar and a table listing news items. The table has columns for ID, Title, Create Date, Update Date, Deleted, and Modify. A single row is visible with the following data: ID 1, Title Test, Create Date 2018-05-17, Update Date 2018-05-17, Deleted N, and Modify (a blue button). The 'Modify' button is highlighted with a red box. On the left side of the table, there's a large empty box labeled 'Contents'. At the bottom of the screen, there are navigation icons and a 'Save' button.

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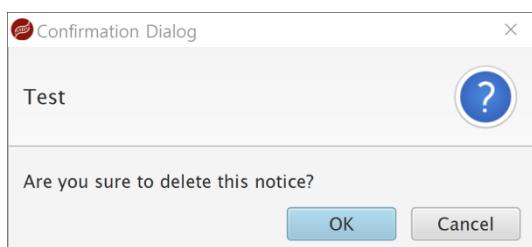
The screenshot shows the 'System Logs' page. At the top, there is a navigation bar with several tabs: Analysis Status, User Account, Panel, Report, News&Tips, **System Logs**, and Clinical Variant DB. The 'System Logs' tab is highlighted with a red box and a yellow circle containing the number 1. Below the navigation bar, there is a search and filter section. It includes fields for Date (2018-05-17 17:27:05), Log Type (Login), Log Msg, User (Administrator), and Group (Administrators). To the right of these fields is a 'Refresh' button (3) enclosed in a red box with a yellow circle. At the bottom of the page, there is a pagination control with buttons for <<, <, 1, >, >>, and a text input for 'Page' followed by a 'GO' button.

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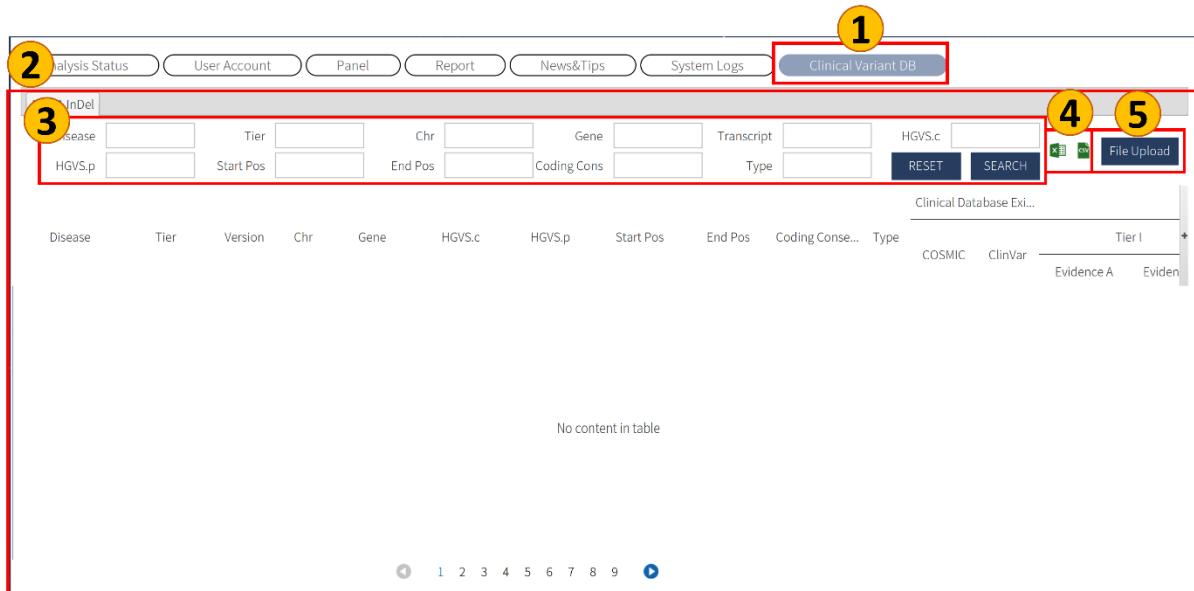The screenshot shows the 'Clinical Variant DB' interface. At the top, there are several tabs: Analysis Status, User Account, Panel, Report, News&Tips, System Logs, and Clinical Variant DB. Below these are search fields for Disease, Tier, Start Pos, End Pos, Chr, Gene, HGVS.c, HGVS.p, Transcript, Coding Cons, Type, HGVS.c, HGVS.p, Start Pos, End Pos, Coding Consequence, Type, COSMIC, ClinVar, and Tier I. A 'File Upload' button is located at the top right of the search area, highlighted with a red box. Below the search area is a table header with columns: Disease, Tier, Version, Chr, Gene, HGVS.c, HGVS.p, Start Pos, End Pos, Coding Consequence, Type, COSMIC, ClinVar, and Tier I. The main content area below the table header displays the message 'No content in table'. At the bottom of the page are navigation icons for back, forward, and search.

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약어	설명
AA	Amino Acid
ACMG	The American College of Medical Genetics and Genomics
AMP	Association of Molecular Pathology
API	Application Programming Interface
ASCO	The American Society of Clinical Oncology
BAM	Binary sequence Alignment Map
CAP	College of American Pathologists
cDNA	Complementary DNA
CNV	Copy Number Variation
CPIC	The Clinical Pharmacogenetics Implementation Consortium
GUI	Graphical User Interface
HGVS	Human Genome Variation Society
HUGO	Human Genome Organization
IGV	Integrative Genomics Viewer by Broad Institute
Indel	Insertion / Deletion
JSON	JavaScript Object Notation
MAF	Min Allele Fraction
NT	NucleoTide
PMID	PubMed ID
QC	Quality Control
RCV	Reference ClinVar
ROI	Region of Interest
SNP	Single Nucleotide Polymorphism
SNV	Single Nucleotide Variants
UI	User Interface
URL	Uniform Resource Locator
UX	User Experience
VAF	Variant Allele Frequency
VCF	Variant Call Format

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약어	정식 명칭	설명
1KGP	1000 Genomes Project	Database of variants found during low-coverage and high-coverage genomic and targeted sequencing from 26 populations. Provides more diversity compared to the Exome Variant Server but also contains lower-quality data, and some cohorts contain related individuals.
ARUP	ARUP-BRCA Mutation Database	Through collaborations with the University of Utah Huntsman Cancer Institute (HCI), and with the WHO International Agency for Research on Cancer (IARC), the University of Utah Department of Pathology and ARUP Laboratories is pleased to host the BRCA mutation database.
Be	BRCA Exchange	Curated data platform that provides information on variants of specific genes: BRCA1 and BRCA2.
BIC	Breast Cancer Information Core	An Open Access On-Line Breast Cancer Mutation Data Base.
ClinVar	ClinVar	Database of assertions about the clinical significance and phenotype relationship of human variations.
COSMIC	Catalogue Of Somatic Mutations In Cancer	Database for exploring the impact of somatic mutations in human cancer.
dbNSFP	Database for Nonsynonymous SNPs' Functional Predictions	Database developed for functional prediction and annotation of all potential non-synonymous single-nucleotide variants (nsSNVs) in the human genome.
dbscSNV	Database for SNVs within splicing consensus regions	dbscSNV includes all potential human SNVs within splicing consensus regions (-3 to +8 at the 5' splice site and -12 to +2 at the 3' splice site), i.e. scSNVs, related functional annotations and two ensemble prediction scores for predicting their potential of altering splicing.
dbSNP	the NCBI database of genetic variation	Database of short genetic variations (typically = <50bp) submitted from many sources. May lack details of the originating study and may contain pathogenic variants.

FATHMM	Functional Analysis through Hidden Markov Models	A high-throughput web-server capable of predicting the functional consequences of both coding variants, i.e. non-synonymous single nucleotide variants (nsSNVs), and non-coding variants.
ENIGMA	Evidence-based Network for the Interpretation of Germline Mutant Alleles	Evidence-based Network for the Interpretation of Germline Mutant Alleles.
ESP6500	NHLBI Grand Opportunity Exome Sequencing Project	Database of variants found during exome sequencing of several large cohorts of individuals of European and African American ancestry. Includes coverage data to inform the absence of variation. (6,503 samples drawn from multiple ESP cohorts and represents all of the ESP exome variant data).
ExAC	Exome Aggregation Consortium	Database of variants found during exome sequencing of 61,486 unrelated individuals sequenced as part of various disease-specific and population genetic studies. Pediatric disease subjects as well as related individuals were excluded.
GERP	Genomic Evolutionary Rate Profiling	GERP identifies constrained elements in multiple alignments by quantifying substitution deficits.
gnomAD	Genome Aggregation Database	The Genome Aggregation Database (gnomAD), is a coalition of investigators seeking to aggregate and harmonize exome and genome sequencing data from a variety of large-scale sequencing projects, and to make summary data available for the wider scientific community. The data set provided on this website spans 123,136 exomes and 15,496 genomes from unrelated individuals sequenced as part of various disease-specific and population genetic studies.
KoEXID	Korean Exome Information Database	Analyzed data obtained by deep sequencing 1,303 Korean whole exomes; the data were generated by three independent whole exome sequencing projects (named the KOEX study).

KOHBRA	The Korean Hereditary Breast Cancer	In 2007, the Korean Hereditary Breast Cancer (KOHBRA) Study was established to obtain evidence for the accurate risk assessment and management of hereditary breast and ovarian cancer (HBOC) in Korea. Breast Cancer Res Treat (2015)
KOVA	Korean Variant Archive	To establish a genetic database of healthy Korean individuals, WES data from normal tissues from cancer patients (675 samples; 472 blood and 203 adjacent normal tissues) and blood samples from healthy individuals with no apparent clinical history (380 samples) were collected.
KRGDB	Korean Reference Genome	Since 2012, Center for Genome Science of Korea National Research Institute of Health has conducted the whole genome sequencing project for 1,722 Korean individuals. The project investigated and analyzed 622 samples until 2014 (the first step), and 1,100 samples until 2016 (the second step).
LRT	Likelihood Ratio Test	A statistical test used for comparing the goodness of fit of two statistical models
MetaSVM	Meta-analytic Support Vector Machine	MetaSVM accommodate multiple omics data, making it possible to detect consensus genes associated with diseases across studies.
Mutation Assessor	Mutation Assessor	A database and scoring system to predict the functional impact of mutations in translated proteins in human cancers
Mutation Taster	Mutation Taster	Mutation prediction for the deep-sequencing age
SIFT	Scale-Invariant Feature Transform	Predicting amino acid changes that affect protein function

PATHOGENICITY CLASS

Pathogenic

- Nonsense mutations (stop_gained) with 1KGP and ESP6500 <= 0.01.
- Frameshift (frameshift_variant) with 1KGP and ESP6500 <= 0.01.
- Splice acceptor (splice_acceptor_variant) sites (-2 and -1) and/or splice donor (splice_donor_variant) sites (+1 and +2).
- Loss of natural initiation of transcription (start_codon, start_lost) or termination (stop_codon, stop_lost).
- In frame INDELs (inframe_insertion, inframe_deletion) of more than 3 nucleotides.
- Mutation reported as pathogenic in BIC or ClinVar.

Likely Pathogenic

- Missense (missense_variant) mutation with 1KGP and ESP6500 < 0.01.
- In frame INDELs of 3 nucleotides.
- Missense or synonymous with at least 2 values of predictive score (SIFT, PolyPhen2, MutationTaster) > 0.9.

Likely Benign

- None.

Benign

- Missense, synonymous or intronic variant with 1KGP or ESP6500 > 0.05 with 2 predictive.

Uncertain Significance

- Mutation no classified in pathogenic, likely pathogenic and benign.

TIER CLASS

Tier I

- Variants with strong clinical significance.
- Level A: biomarkers that predict response or resistance to US FDA-approved therapies for a specific type of tumor or have been included in professional guidelines as therapeutic, diagnostic, and/or prognostic biomarkers for specific types of tumors.
- Level B: biomarkers that predict response or resistance to a therapy based on well-powered studies with consensus from experts in the field, or have diagnostic and/or prognostic significance of certain diseases based on well powered studies with expert consensus.
- Variants with strong clinical significance.

Tier II

- Variants with potential clinical significance.
- Level C: biomarkers that predict response or resistance to therapies approved by FDA or professional societies for a different tumor type (ie, off-label use of a drug), serve as inclusion criteria for clinical trials, or have diagnostic and/or prognostic significance based on the results of multiple small studies.
- Level D: biomarkers that show plausible therapeutic significance based on preclinical studies, or may assist disease diagnosis and/or prognosis themselves or along with other biomarkers based on small studies or multiple case reports with no consensus.

Tier III

- Variants of unknown clinical significance.
- Not observed at a significant allele frequency in the general or specific subpopulation databases, or pan-cancer or tumor-specific variant databases.
- No convincing published evidence of cancer association.

Tier IV

- Variants deemed benign or likely benign.
- Observed at significant allele frequency in the general or specific subpopulation databases.
- No existing published evidence of cancer association.

PATHOGENICITY CLASS - ClinVar

Clinical significance on ClinVar submissions (SCV)

Clinical Significance Value	Guidance for Use in ClinVar SCV Records
Benign	As recommended by ACMG/AMP for variants interpreted for Mendelian disorders.
Likely benign	As recommended by ACMG/AMP for variants interpreted for Mendelian disorders.
Uncertain significance	As recommended by ACMG/AMP for variants interpreted for Mendelian disorders.
Likely pathogenic	As recommended by ACMG/AMP for variants interpreted for Mendelian disorders.
Pathogenic	As recommended by ACMG/AMP for variants interpreted for Mendelian disorders.
drug response	A general term for a variant that affects a drug response, not a disease. We anticipate adding more specific drug response terms based on a recommendation by CPIC.
association	For variants identified in a GWAS study and further interpreted for their clinical significance.
risk factor	For variants that are interpreted not to cause a disorder but to increase the risk.
protective	For variants that decrease the risk of a disorder, including infections.
Affects	For variants that cause a non-disease phenotype, such as lactose intolerance.
conflicting data from submitters	Only for submissions from a consortium, where groups within the consortium have conflicting interpretations of a variant but provide a single submission to ClinVar.
other	If ClinVar does not have the appropriate term for your submission, we ask that you submit "other" as clinical significance and contact us to discuss if there are other terms we should add.
not provided	For submissions without an interpretation of clinical significance. The primary goal of ClinVar is to archive reports of clinical significance of variants. Therefore, submissions with a clinical significance of "not provided" should be limited to: <ul style="list-style-type: none"> • "literature only" submissions that report a publication about the variant, without interpreting the clinical significance. • "research" submissions that provide functional significance (e.g.

	<p>undetectable protein level) but no interpretation of clinical significance.</p> <ul style="list-style-type: none"> “clinical testing” or “phenotyping only” submissions from clinics or physicians that provide additional information about individuals with the variant, such as observed phenotypes, but do not interpret the clinical significance.
--	---

(<https://www.ncbi.nlm.nih.gov/clinvar/docs/clinsig/>)

Clinical significance on ClinVar aggregate records (RCV)

Combination of Values from Different Submitters	Reported as
Pathogenic and Likely pathogenic	Pathogenic/Likely Pathogenic
(Pathogenic or Likely pathogenic or Benign or Likely benign) AND Uncertain significance	Conflicting interpretations of pathogenicity
(Pathogenic or Likely pathogenic) AND (Benign or Likely benign)	Conflicting interpretations of pathogenicity
Benign and Likely benign	Benign/Likely benign
Any single ACMG value AND any non-ACMG value, e.g. Uncertain significance and risk factor	Uncertain significance, risk factor
Conflicting ACMG values AND any non-ACMG value, e.g. Pathogenic and Uncertain significance and risk factor	Conflicting interpretations of pathogenicity, risk factor
No ACMG value AND multiple non-ACMG values, e.g. Drug response and risk factor	drug response, risk factor
A consortium does not agree on which ACMG assessment is valid	conflicting data from submitters

(<https://www.ncbi.nlm.nih.gov/clinvar/docs/clinsig/>)

PATHOGENICITY CLASS - BIC

Class	Clinical Significance Value	Description
Class 1	Not pathogenic low clinical significance	There is significant evidence against this variant being a dominant high-risk pathogenic mutation.
Class 2	Likely Notice pathogenic little clinical significance	There is evidence against this variant being a dominant high-risk pathogenic mutation.
Class 3	Uncertain	There is insufficient evidence to place this variant in Class 1,2,4 or 5.
Class 4	Likely Pathogenic	There is evidence that this variant is a dominant high-risk pathogenic mutation.
Class 5	Pathogenic	There is significant evidence to suggest that this variant is a dominant high-risk pathogenic mutation.

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필드 명	설명
사용자 설정 항목	
Pathogenicity	Pathogenic class using ACMG guideline.
Warning	Warning variant.
Report	Reported variant.
변이 및 유전자 정보	
Gene	Gene Symbol (HUGO Gene Nomenclature).
Strand	Whether the gene is located on the plus or minus strand.
Transcript Accession	NCBI Transcript Accession ID.
Type	Variant type SNP: Single-Nucleotide Polymorphism, INS: Insertion, Del: Deletion, MNP: Multiple-nucleotide polymorphism, Complex: Multiple-nucleotide and an Indel.
Consequence	Sequence change result by variant.
NT Change	The variant coordinates relative to the coding DNA
NT Change (BIC)	The variant coordinates relative to the coding DNA (BRCA format- BRCA1: A of start codon ATG is #120, BRCA2: A of start codon ATG is #229).
Protein Accession	NCBI Protein Accession ID.
AA Change	Change amino acid. (three letter)
AA Change (Single)	Change amino acid (single letter).
변이 위치 정보	
Chr	The chromosome numbers.
Start Position	The variant coordinate on the reference genome.
Ref	Reference allele.
Alt	Alternative allele.
Allele 정보	
Fraction	The proportion of the reads supporting the variant count.
Depth	Total coverage of the variant.
Ref Count	Number of reads supporting the reference allele.
Alt Count	Number of reads supporting the alternative allele.
Exon	Number of exon in which the variant is located exon/total exon.

Exon (BIC)	Exon number; Reference sequences: BRCA1 GenBank U14680; BRCA2 GenBank U43746.
dbSNP ID	dbSNP RS ID.
dbSNP Common ID	dbSNP RS Common ID.
Zygosity	The words homozygous, heterozygous, and hemizygous are used to describe the genotype of a diploid organism at a single locus on the DNA.
ClinVar 데이터베이스 정보	
ClinVar Variation ID	ClinVar variation ID
ClinVar Interpretation	ClinVar variation interpretation
ClinVar Review Status	ClinVar review status for the Variation ID
ClinVar Submitted Accession	ClinVar RCV accessions.
ClinVar Submitted Interpretation	Clinical classification.
ClinVar Submitted Disease	ClinVar disease name.
ClinVar Submitted OMIM	OMIM number reported to ClinVar
In silico 데이터베이스 정보	
SIFT Prediction	SIFT's prediction value
Mutation Taster Prediction	Mutation Taster's mutation prediction
Mutation Assessor Prediction	Mutation Assessor's driver mutation prediction
LRT Prediction	LRT's prediction value
FATHMM Prediction	FATHMM's the functional consequences prediction value
MetaSVM Prediction	Prediction of our SVM based ensemble prediction score
dbscSNV ADA Score	Adaboost classifier score from dbscSNV
dbscSNV RF Score	Random forest classifier score from dbscSNV
GERP++ NR	GERP's the corresponding neutral rate
GERP++ RS	GERP's rejected substitutions
Population 빈도	
1KGP ALL	1KGP's All population frequency.
1KGP AFR	1KGP's African population frequency.
1KGP AMR	1KGP's American population frequency.
1KGP EAS	1KGP's East Asian population frequency.
1KGP EUR	1KGP's European population frequency.
1KGP SAS	1KGP's South Asian population frequency.

ESP All	ESP6500's All population frequency.
ESP AA	ESP6500's African American population frequency.
ESP EA	ESP6500's European American population frequency.
ExAC	ExAC population frequency.
KRGDB	KRGDB population frequency.
KoEXID	KoEXID population frequency.
BRCA 데이터베이스 정보	
BIC Category	Mutation type in BIC. 3UTR: 3'UTR. 5TR: 5' UTR. F: Frameshift. IFI: In frame insertion. IFD: In frame deletion. IVS: Intervening sequence. M: Missense. N: Nonsense. P: Polymorphism. S: Splice. Syn: Synonymous. UV: Unclassified variant.
BIC Class	Clinical classification using 5 classes in BIC. 0: pending. 1: Not pathogenic/low clinical significance. 2: Likely not pathogenic/little clinical significance. 3: Uncertain. 4: Likely pathogenic. 5: Pathogenic.
BIC Designation	Designation of described mutation according to BIC nomenclature guidelines.
BIC Importance	Clinical importance in BIC.
BIC NT	Nucleotide # in the transcript (cDNA) at which mutation occurs; Reference sequences: BRCA1 GenBank U14680; BRCA2 GenBank U43746.
KOHBRA Frequency	High risk breast cancer frequency in KOHBRA.
KOHBRA Patient	High risk breast cancer patient in KOHBRA.
Be BIC Category	Mutation type in BIC.
Be BIC Ethnic	Ethnic of BIC.
Be BIC Nationality	Nationality of the patient associated with the variant. From BIC.
Be BIC Pathogenicity	BIC pathogenicity.

Be ClinVar Method	Method used to collect the data that supports the assertion of clinical significance. Allowed values: case-control, clinical testing, literature only, reference population, research.
Be ClinVar Origin	Variants are classified as either germline or somatic, depending on how they are acquired. Germline variants are genetic changes that we inherit from our parents. Somatic variants are DNA changes that we acquire over our lifetime, often through exposure to pollutants, toxins, radiation and other carcinogens.
Be ClinVar Pathogenicity	ClinVar pathogenicity.
Be ClinVar Update	Date the variant was last updated in ClinVar.
Be ENIGMA Condition	Human-readable condition, describing the biological impact of the variation.
Be ENIGMA Pathogenicity	IRAC Classification from exLOVD.
Be ENIGMA Update	The date on which the clinical significance of the variant was last evaluated by ENIGMA.
Be Gene	Reference gene.
Be NT	Nucleotide change.
Be Transcript	Reference transcript.
ENIGMA	The Clinical Significance columns indicates whether expert curators have determined if the variant is pathogenic or benign. What do these classifications mean? Pathogenic: variants confer an increased risk of disease. Likely pathogenic: variants have good evidence to support an association with disease risk. Likely benign: variants have good evidence to support no association with disease risk. Benign: variants are not associated with any markedly increased risk of disease. Variants of uncertain significance (VUS): are those for which the evidence of disease risk is not clear yet, sometimes because there is not yet enough evidence to classify them as either pathogenic or benign.

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필드 명	설명
사용자 설정 항목	
Tier	Tier class using AMP/ASCO/CAP guideline.
Warning	Warning variant.
FALSE	False positive variant.
Report	Reported variant.
변이 및 유전자 정보	
Gene	Gene Symbol (HUGO Gene Nomenclature).
Strand	Whether the gene is located on the plus or minus strand.
Transcript Accession	NCBI Transcript Accession ID.
Type	Variant type SNP: Single-Nucleotide Polymorphism, INS: Insertion, Del: Deletion, MNP: Multiple-nucleotide polymorphism, Complex: Multiple-nucleotide and an Indel.
Type Extension	Variant clinical type
Consequence	Sequence change result by variant.
NT Change	The variant coordinates relative to the coding DNA
Protein Accession	NCBI Protein Accession ID.
AA Change	Change amino acid. (three letter)
AA Change (Single)	Change amino acid (single letter).
변이 위치 정보	
Chr	The chromosome numbers.
Start Position	The variant coordinate on the reference genome.
Ref	Reference allele.
Alt	Alternative allele.
Allele 정보	
Fraction	The proportion of the reads supporting the variant count.
Depth	Total coverage of the variant.
Ref Count	Number of reads supporting the reference allele.
Alt Count	Number of reads supporting the alternative allele.
Exon	Number of exon in which the variant is located exon/total exon.
dbSNP ID	dbSNP RS ID.
dbSNP Common ID	dbSNP RS Common ID.
COSMIC ID	COSMIC ID.
COSMIC Occurrence	Tissue classification in COSMIC.
COSMIC Count	Sample count detected COSMIC variant.
ClinVar 데이터베이스 정보	

ClinVar Variation ID	ClinVar variation ID
ClinVar Interpretation	ClinVar variation interpretation
ClinVar Review Status	ClinVar review status for the Variation ID
ClinVar Submitted Accession	ClinVar RCV accessions.
ClinVar Submitted Interpretation	Clinical classification.
ClinVar Submitted Disease	ClinVar disease name.
ClinVar Submitted OMIM	OMIM number reported to ClinVar
<i>In silico</i> 데이터베이스 정보	
SIFT Prediction	SIFT's prediction value
Mutation Taster Prediction	Mutation Taster's mutation prediction
Mutation Assessor Prediction	Mutation Assessor's driver mutation prediction
LRT Prediction	LRT's prediction value
FATHMM Prediction	FATHMM's the functional consequences prediction value
GERP++ NR	GERP's the corresponding neutral rate
GERP++ RS	GERP's rejected substitutions
Population 빈도	
1KGP ALL	1KGP's All population frequency.
1KGP AFR	1KGP's African population frequency.
1KGP AMR	1KGP's American population frequency.
1KGP EAS	1KGP's East Asian population frequency.
1KGP EUR	1KGP's European population frequency.
1KGP SAS	1KGP's South Asian population frequency.
ESP All	ESP6500's All population frequency.
ESP AA	ESP6500's African American population frequency.
ESP EA	ESP6500's European American population frequency.
ExAC	ExAC population frequency.
gnomAD All	gnomAD's All population frequency.
gnomAD AFR	gnomAD's African American population frequency.
gnomAD AMR	gnomAD's Admixed American population frequency.
gnomAD EAS	gnomAD's East Asian population frequency.
gnomAD FIN	gnomAD's Finnish population frequency.
gnomAD NFE	gnomAD's Non Finnish European population frequency.
gnomAD OTH	gnomAD's Others population frequency.

gnomAD SAS	gnomAD's South Asian population frequency.
KRGDB	KRGDB population frequency.
KoEXID	KoEXID population frequency.

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필드 명	설명
사용자 설정 항목	
Tier	Tier class using AMP/ASCO/CAP guideline.
Warning	Warning variant.
FALSE	False positive variant.
Report	Reported variant.
변이 및 유전자 정보	
Gene	Gene Symbol (HUGO Gene Nomenclature).
Strand	Whether the gene is located on the plus or minus strand.
Transcript Accession	NCBI Transcript Accession ID.
Type	Variant type SNP: Single-Nucleotide Polymorphism, INS: Insertion, Del: Deletion, MNP: Multiple-nucleotide polymorphism, Complex: Multiple-nucleotide and an Indel.
Type Extension	Variant clinical type
Consequence	Sequence change result by variant.
NT Change	The variant coordinates relative to the coding DNA
Protein Accession	NCBI Protein Accession ID.
AA Change	Change amino acid. (three letter)
AA Change (Single)	Change amino acid (single letter).
변이 위치 정보	
Chr	The chromosome numbers.
Start Position	The variant coordinate on the reference genome.
Ref	Reference allele.
Alt	Alternative allele.
Allele 정보	
Fraction	The proportion of the reads supporting the variant count.
Depth	Total coverage of the variant.
Ref Count	Number of reads supporting the reference allele.
Alt Count	Number of reads supporting the alternative allele.
Exon	Number of exon in which the variant is located exon/total exon.
dbSNP ID	dbSNP RS ID.

dbSNP Common ID	dbSNP RS Common ID.
COSMIC ID	COSMIC ID.
COSMIC Occurrence	Tissue classification in COSMIC.
COSMIC Count	Sample count detected COSMIC variant.
ClinVar 데이터베이스 정보	
ClinVar Variation ID	ClinVar variation ID
ClinVar Interpretation	ClinVar variation interpretation
ClinVar Review Status	ClinVar review status for the Variation ID
ClinVar Submitted Accession	ClinVar RCV accessions.
ClinVar Submitted Interpretation	Clinical classification.
ClinVar Submitted Disease	ClinVar disease name.
ClinVar Submitted OMIM	OMIM number reported to ClinVar
Population 빈도	
1KGP ALL	1KGP's All population frequency.
1KGP AFR	1KGP's African population frequency.
1KGP AMR	1KGP's American population frequency.
1KGP EAS	1KGP's East Asian population frequency.
1KGP EUR	1KGP's European population frequency.
1KGP SAS	1KGP's South Asian population frequency.
ESP All	ESP6500's All population frequency.
ESP AA	ESP6500's African American population frequency.
ESP EA	ESP6500's European American population frequency.
ExAC	ExAC population frequency.
gnomAD All	gnomAD's All population frequency.
gnomAD AFR	gnomAD's African American population frequency.
gnomAD AMR	gnomAD's Admixed American population frequency.
gnomAD EAS	gnomAD's East Asian population frequency.
gnomAD FIN	gnomAD's Finnish population frequency.
gnomAD NFE	gnomAD's Non Finnish European population frequency.
gnomAD OTH	gnomAD's Others population frequency.
gnomAD SAS	gnomAD's South Asian population frequency.
KRGDB	KRGDB population frequency.
KoEXID	KoEXID population frequency.

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필드 명	설명
사용자 설정 항목	
Pathogenicity	Pathogenic class using ACMG guideline.
Warning	Warning variant.
Common Variants	Common SNV/InDel in Population DB or in House DB
FALSE	False positive variant.
Report	Reported variant.
변이 및 유전자 정보	
Gene	Gene Symbol (HUGO Gene Nomenclature).
Inheritance	Gene Inheritance
Strand	Whether the gene is located on the plus or minus strand.
Transcript Accession	NCBI Transcript Accession ID.
Type	Variant type SNP: Single-Nucleotide Polymorphism, INS: Insertion, Del: Deletion, MNP: Multiple-nucleotide polymorphism, Complex: Multiple-nucleotide and an Indel.
Type Extension	Variant clinical type
Consequence	Sequence change result by variant.
NT Change	The variant coordinates relative to the coding DNA
Protein Accession	NCBI Protein Accession ID.
AA Change	Change amino acid. (three letter)
AA Change (Single)	Change amino acid (single letter).
변이 위치 정보	
Chr	The chromosome numbers.
Start Position	The variant coordinate on the reference genome.
Ref	Reference allele.
Alt	Alternative allele.
Allele 정보	
Fraction	The proportion of the reads supporting the variant count.
Depth	Total coverage of the variant.
Ref Count	Number of reads supporting the reference allele.
Alt Count	Number of reads supporting the alternative allele.
Exon	Number of exon in which the variant is located exon/total exon.
dbSNP ID	dbSNP RS ID.
Zygosity	The words homozygous, heterozygous, and hemizygous are used to describe the genotype of a diploid organism at a single locus on the DNA.

ClinVar 데이터베이스 정보	
ClinVar Variation ID	ClinVar variation ID
ClinVar Interpretation	ClinVar variation interpretation
ClinVar Review Status	ClinVar review status for the Variation ID
ClinVar Submitted Accession	ClinVar RCV accessions.
ClinVar Submitted Interpretation	Clinical classification.
ClinVar Submitted Disease	ClinVar disease name.
ClinVar Submitted OMIM	OMIM number reported to ClinVar
In silico 데이터베이스 정보	
SIFT Prediction	SIFT's prediction value
Mutation Taster Prediction	Mutation Taster's mutation prediction
Mutation Assessor Prediction	Mutation Assessor's driver mutation prediction
LRT Prediction	LRT's prediction value
FATHMM Prediction	FATHMM's the functional consequences prediction value
MetaSVM Prediction	Prediction of our SVM based ensemble prediction score
dbscSNV ADA Score	Adaboost classifier score from dbscSNV
dbscSNV RF Score	Random forest classifier score from dbscSNV
GERP++ NR	GERP's the corresponding neutral rate
GERP++ RS	GERP's rejected substitutions
Population 빈도	
1KGP ALL	1KGP's All population frequency.
1KGP AFR	1KGP's African population frequency.
1KGP AMR	1KGP's American population frequency.
1KGP EAS	1KGP's East Asian population frequency.
1KGP EUR	1KGP's European population frequency.
1KGP SAS	1KGP's South Asian population frequency.
ESP All	ESP6500's All population frequency.
ESP AA	ESP6500's African American population frequency.
ESP EA	ESP6500's European American population frequency.
ExAC	ExAC population frequency.
gnomAD All	gnomAD's All population frequency.
gnomAD AFR	gnomAD's African American population frequency.
gnomAD AMR	gnomAD's Admixed American population frequency.

gnomAD EAS	gnomAD's East Asian population frequency.
gnomAD FIN	gnomAD's Finnish population frequency.
gnomAD NFE	gnomAD's Non Finnish European population frequency.
gnomAD OTH	gnomAD's Others population frequency.
gnomAD SAS	gnomAD's South Asian population frequency.
KRGDB	KRGDB population frequency.
KoEXID	KoEXID population frequency.

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용어	설명
fragment_length	abs(ref - alt) median fragment length
mapping_error	Variation resulting from regions with similar sequences other than genes (ex. Pseudogene)
mapping_quality	ref - alt median mapping quality
orientation_bias	Orientation bias (in one of the specified artifact mode(s) or complement) seen in one or more samples.
sequencing_error	Variation observed simultaneously with low VAF values(< 15%) in multiple samples (eg, a variation that only appears in a specific direction reading)
snp_candidate	Variation observed in VAF close to homo / hetero (near 50% or near 100%) at the same time in multiple samples but not SNP Population DB information
strand_artifact	Evidence for alt allele comes from one read direction only
t_lod	Tumor does not meet likelihood threshold

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용어	설명
homopolymer	InDel mutation in sequences with one base in succession
repeat_sequence	InDel mutation occurring in a sequence in which a pattern consisting of two or more bases is repeated
sequencing_error	Variation observed simultaneously with low VAF values(< 30%) in multiple samples (eg, a variation that only appears in a specific direction reading)
mapping_error	Variation resulting from regions with similar sequences other than genes (ex. Pseudogene)

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- 1)Lynch HT, Snyder CL, Lynch JF, Riley BD, Rubinstein WS. Hereditary breast-ovarian cancer at the bedside: role of the medical oncologist. *Journal of clinical oncology*. 2003;21(4):740-53.
- 2)Robson M, Offit K. Clinical practice. Management of an inherited predisposition to breast cancer. *The New England journal of medicine*. 2007;357(2):154-62.
- 3)Pal T, Permuth-Wey J, Betts JA, et al. BRCA1 and BRCA2 mutations account for a large proportion of ovarian carcinoma cases. *Cancer*. 2005;104(12):2807-16.
- 4)Kang E, Seong MW, Park SK, et al. The prevalence and spectrum of BRCA1 and BRCA2 mutations in Korean population: recent update of the Korean Hereditary Breast Cancer (KOHBRA) study. *Breast Cancer Res Treat*. 2015;151:157-68.
- 5)Moyer VA. Risk assessment, genetic counseling, and genetic testing for BRCA-related cancer in women: U.S. Preventive Services Task Force recommendation statement. *Annals of internal medicine*. 2014;160(4):271-81.
- 6)Gadzicki D, Evans DG, Harris H, et al. Genetic testing for familial/hereditary breast cancer-comparison of guidelines and recommendations from the UK, France, the Netherlands and Germany. *Journal of community genetics*. 2011;2(2):53-69.
- 7)American Society of Clinical Oncology policy statement update: genetic testing for cancer susceptibility. *Journal of clinical oncology*. 2003;21(12):2397-406.
- 8)Rechards S, Aziz N, Bale S, et al. Standards and guidelines for the interpretation of sequence variants: a joint consensus recommendation of the American College of Medical Genetics and Genomics and the Association for Molecular Pathology. *Genetics in Medicine*. 2015;17(5):405-24.
- 9)CPIC Term Standardization for Clinical Pharmacogenetic Test Results Project.
<https://www.pharmgkb.org/page/cpicTermProject>
- 10)Representation of clinical significance in ClinVar and other variation resources at NCBI.
<https://www.ncbi.nlm.nih.gov/clinvar/docs/clinsig/>

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Copyright © 2016 by NGeneBio Co., Ltd. All Rights Reserved.